

The intense, pulsed positron source EPOS at the Research Centre Dresden-Rossendorf

The EPOS Team and R. Krause-Rehberg


Martin-Luther University, Halle-Wittenberg, Dept. of Physics, 06099 Halle / Germany

Institute of Radiation Physics, Research Centre Dresden-Rossendorf, PO Box 510119, 01314 Dresden / Germany

- New Concept of EPOS
- Progress of the Mono-energetic Positron Beam (MePS)
- Gamma-induced Positron Spectroscopy (GiPS)




Ground plan of the ELBE hall

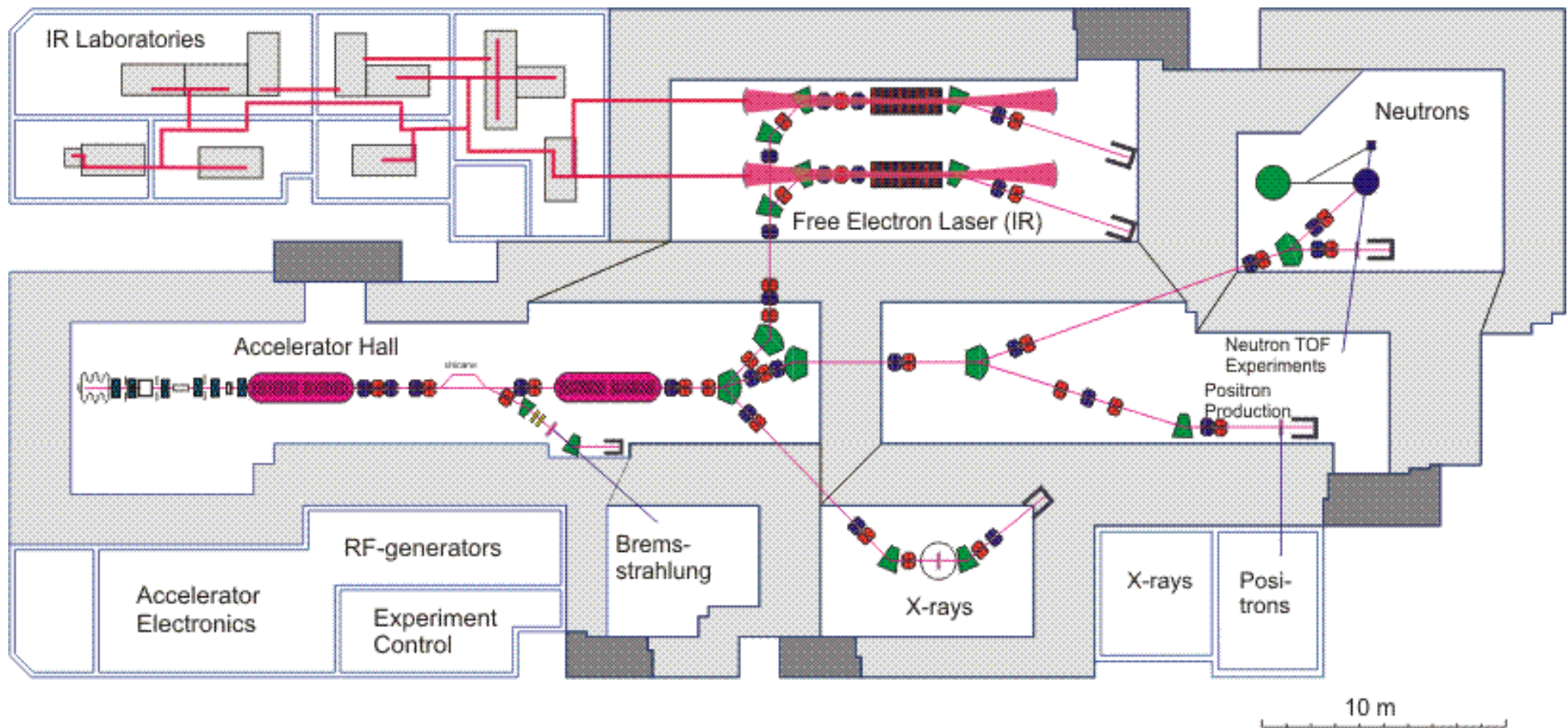

EPOS (ELBE POSitron Source) at ELBE

- Radiation source ELBE = **E**lectron **L**inac with high **B**rilliance and low **E**mittance
- Primary electron beam (40 MeV x 1 mA = 40 kW)
- Time structure: infinite sequence of very short electron bunches (cw-mode)
- Spot size: 5mm diameter

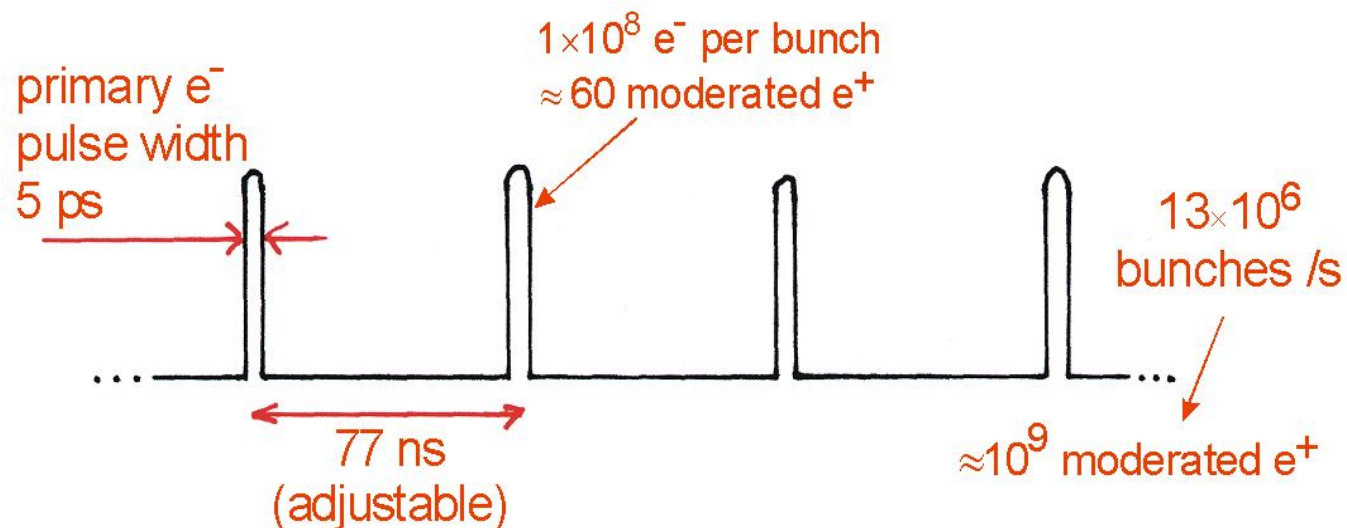

electron bunch structure

New Concept of EPOS (ELBE Positron Source)

MePS

Monoenergetic Positron Spectroscopy

- Cave 111b / Lab 111d
- monoenergetic (slow) positrons
- pulsed system
- LT, CDBS, AMOC
- Still under construction

Information Depth:
0...5 μm

CoPS

Conventional Positron Spectroscopy

- LT, CDBS, AMOC
- using ^{22}Na foil sources
- He-cryostat
- automated system
- digital detector system

Information Depth:
10...200 μm

GiPS


Gamma-induced Positron Spectroscopy

- Cave 109 (nuclear physics)
- Positron generation by Bremsstrahlung
- Information in complete bulky sample (up to 50 cm^3)
- all relevant positron techniques (LT, CDBS, AMOC)

Information Depth:
0.1 mm ...5 cm


Progress of Mono-energetic Positron Beam

- 40 MeV, 1 mA, 26 MHz repetition time in cw mode; lifetime, CDBS and AMOC with slow e^+
- Retain original time structure for simplicity and best time resolution

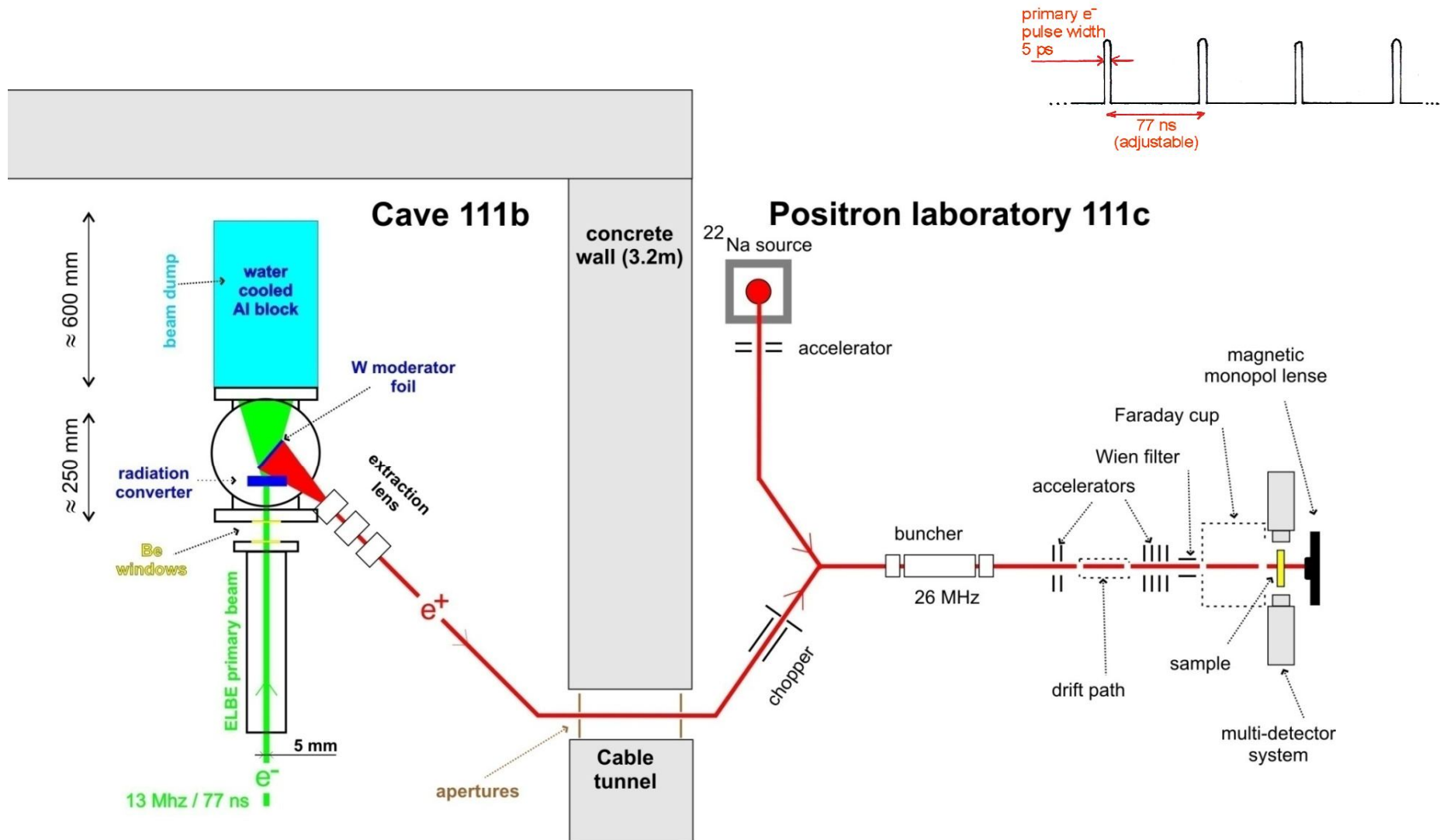


Simulation of Buncher Voltages

Both buncher RF-voltage amplitudes and the drift path energy must be adjusted for each beam energy for optimum time resolution

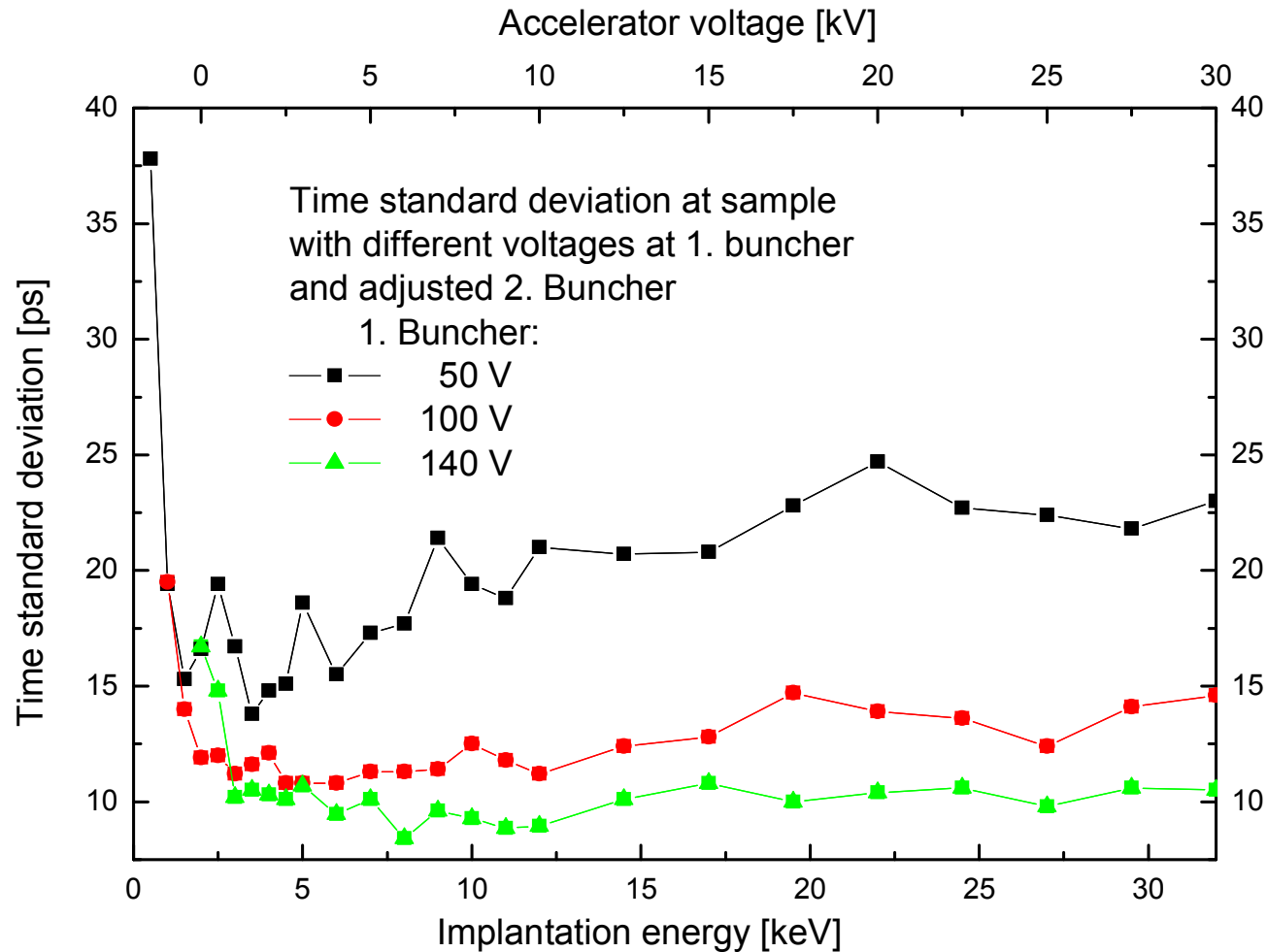


Detector system

- **3 experiments:** lifetime spectroscopy (8 BaF₂ detectors); Doppler coincidence (2 Ge detectors), and AMOC (1 Ge and 1 BaF₂ detector)
- **digital detection system:**
 - lifetime: almost nothing to adjust; time scale exactly the same for all detectors; easy realization of coincidence
 - Doppler: better energy resolution and pile-up rejection expected
 - pulse-shape discrimination improves spectra quality


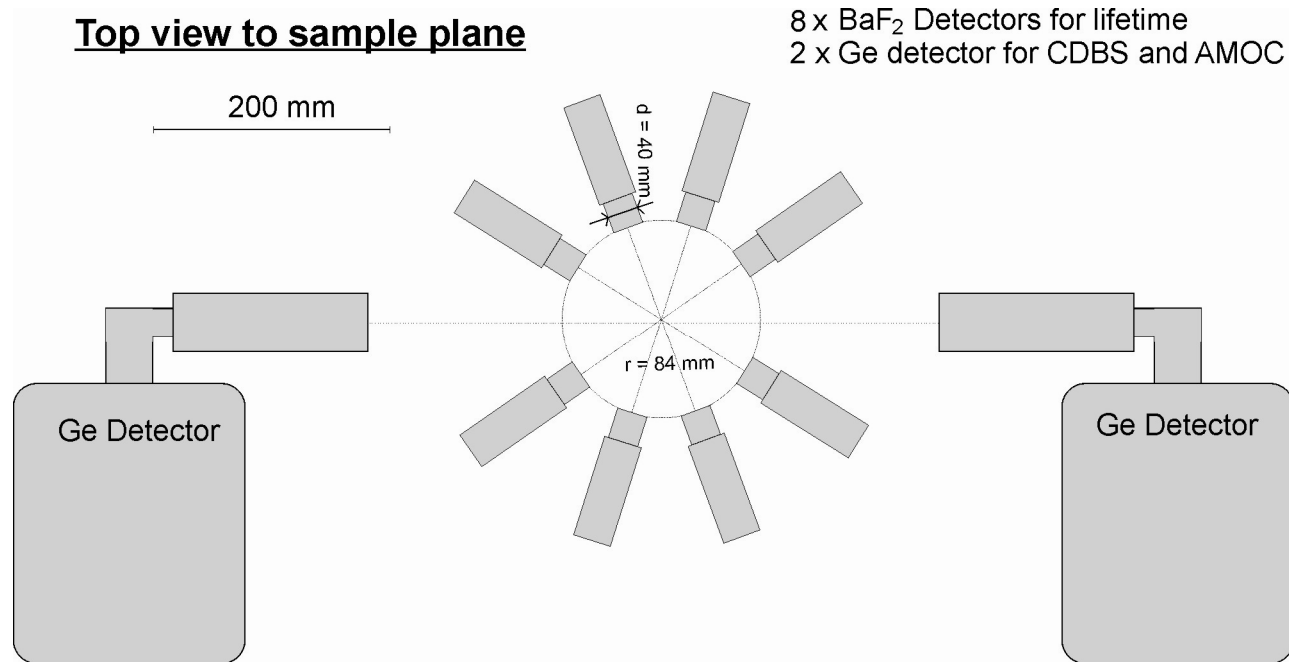


Digital lifetime measurement for ^{22}Na Setup

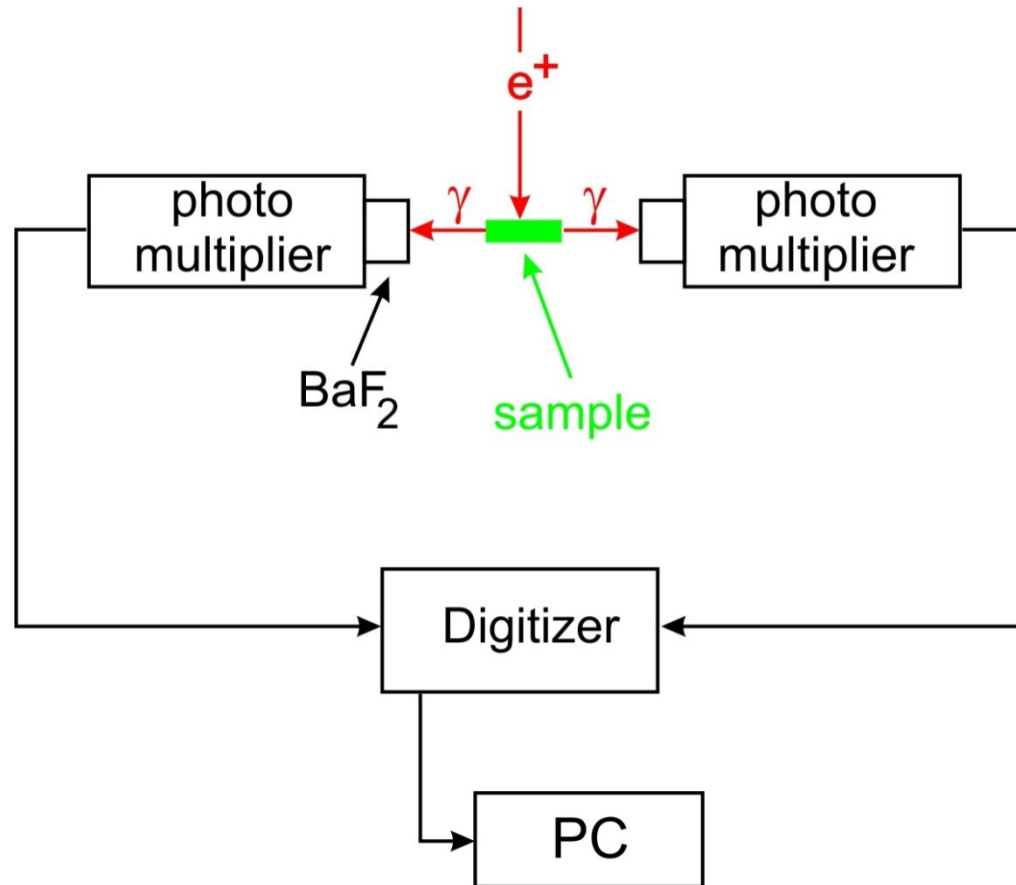


- much simpler setup
- timing very accurate
- pulse-shape discrimination (suppress "bad pulses")
- each detector for start & stop (double statistics)

screenshot of two digitized anode pulses

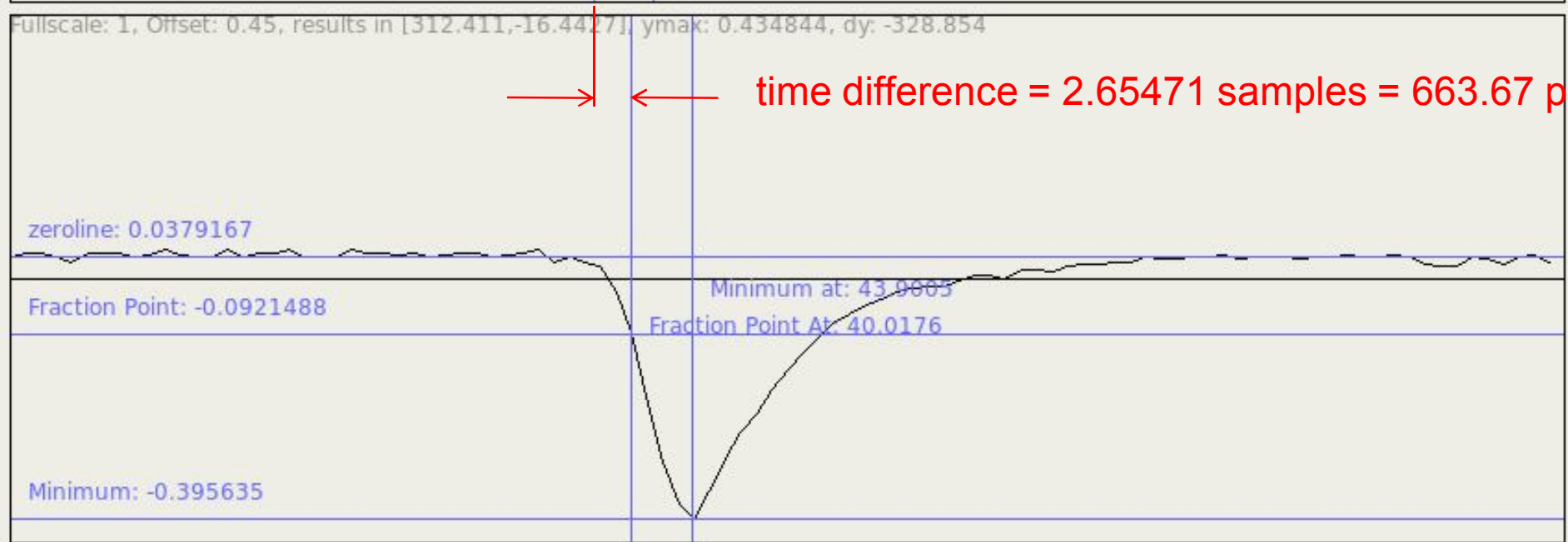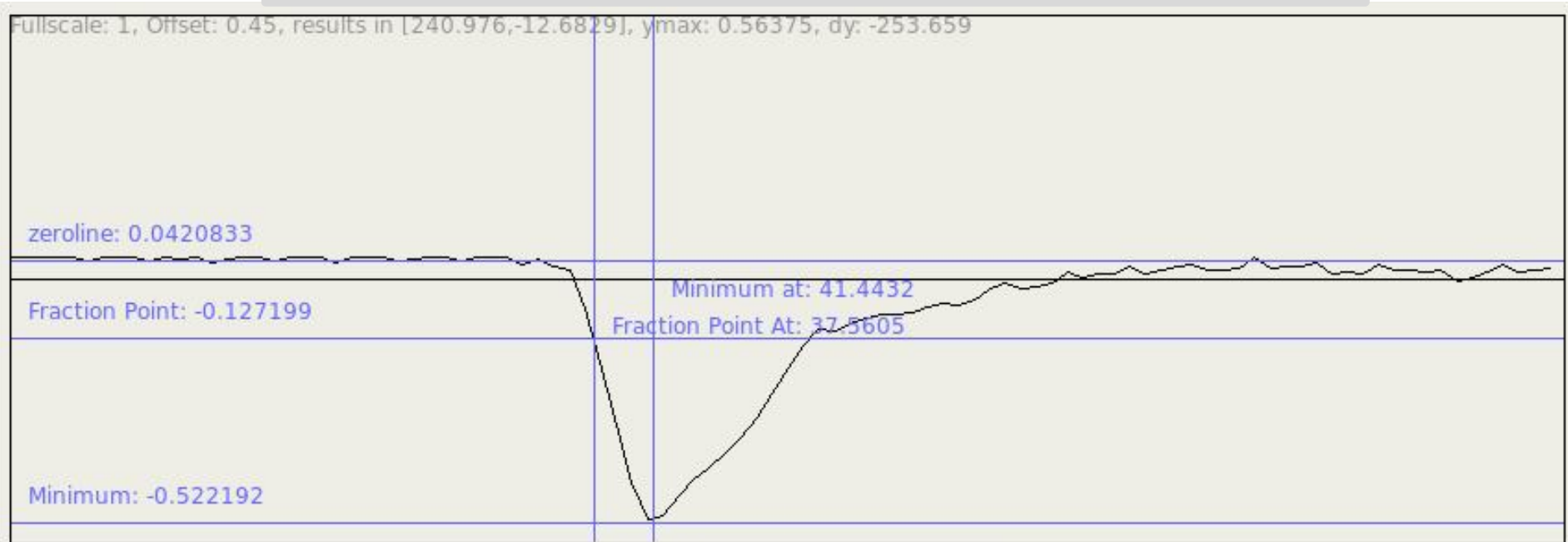

Digital lifetime measurement for EPOS

- coincident registration of both annihilation gamma will reduce background
- will improve time resolution by factor 0.707
- individual registration (no coincidence) will improve count rate

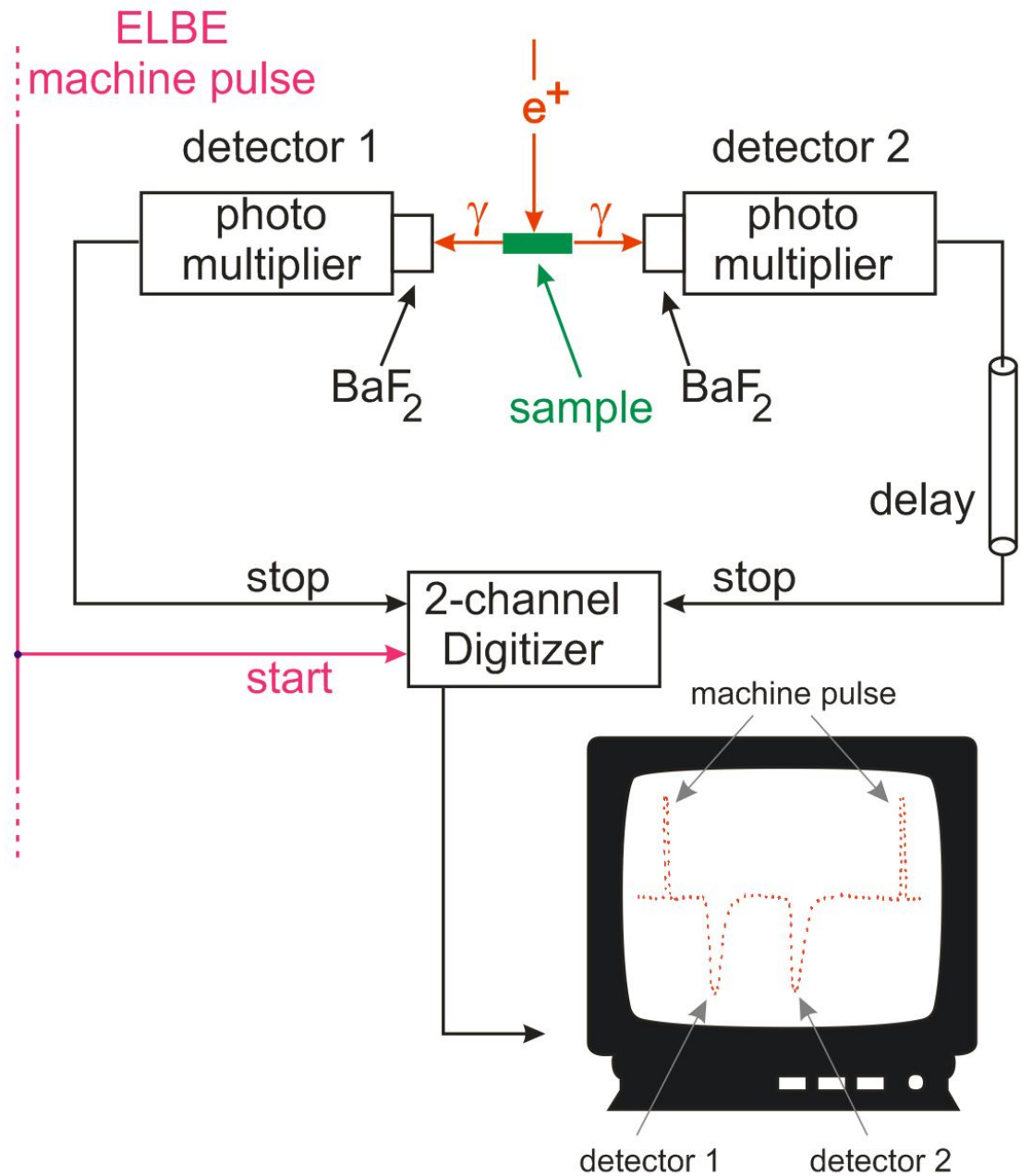



Research Center Dresden-Rossendorf, 28.-30. September 2009

Workshop on Digital Signal Processing in Nuclear Science


<http://positron.physik.uni-halle.de/EPOS/>

Open-source Project


<http://positron.physik.uni-halle.de/EPOS/Software/>

Electron-Positron Converter almost finished



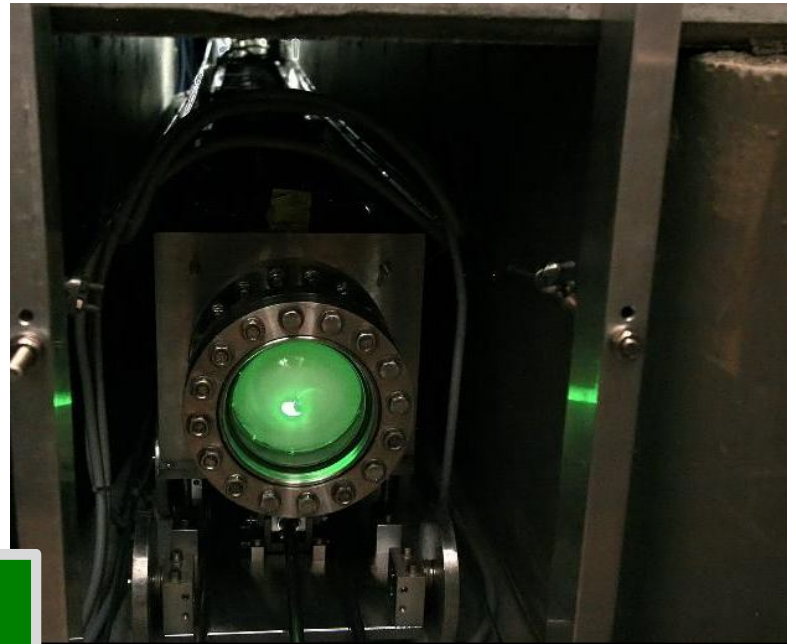

Test of Beam Guidance with an Electron Source

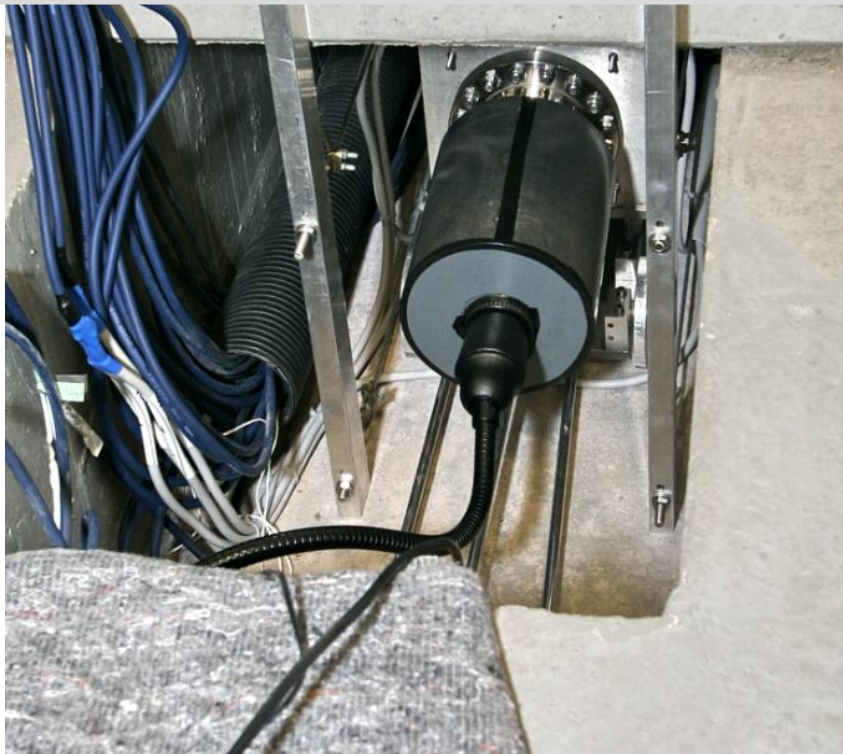


Adjustment of magnetic Guidance System with ^{22}Na source and Moderator





mounting of ^{22}Na source



Still waiting for γ Quanta in Lab ...

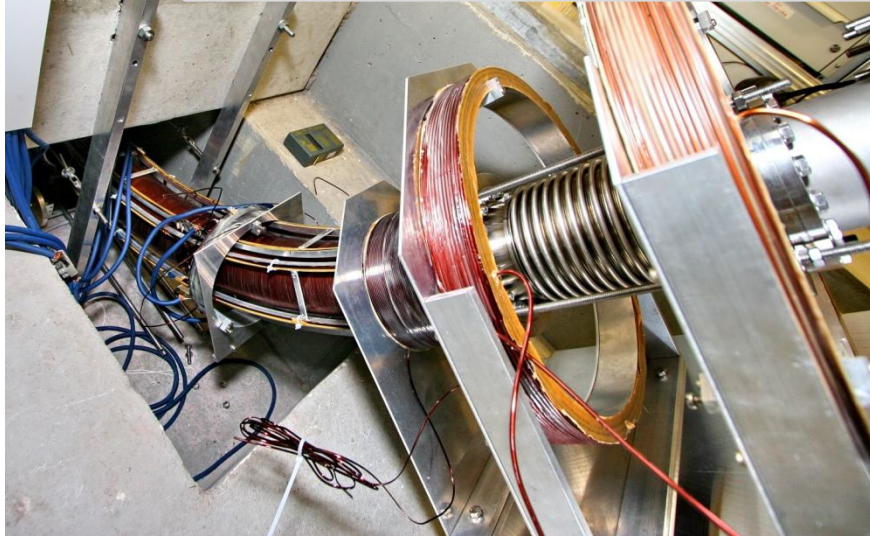

- Problem: 10 x 2 steering coils must be adjusted
- automatic LabView program is looking for annihilation gamma at end of beam line

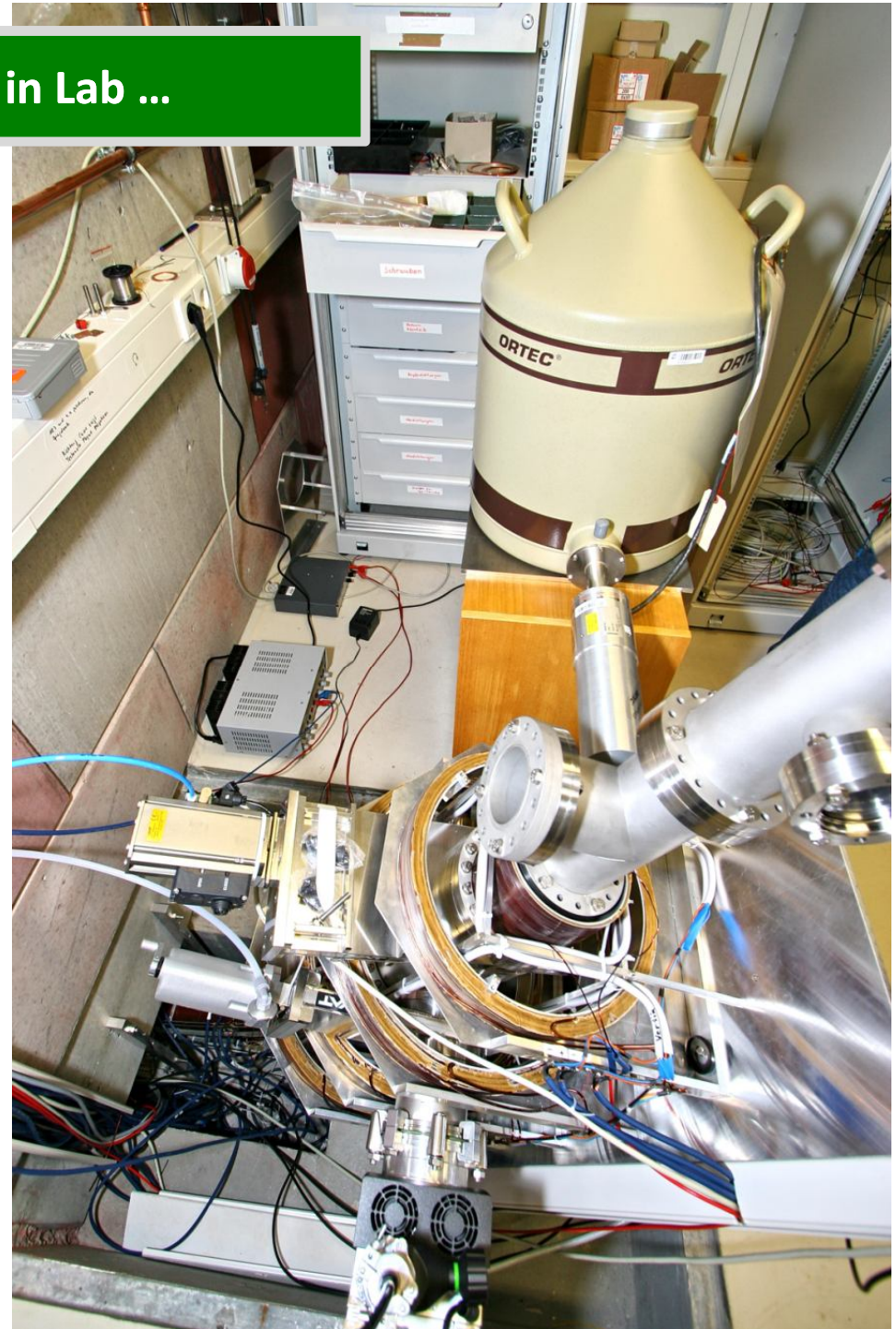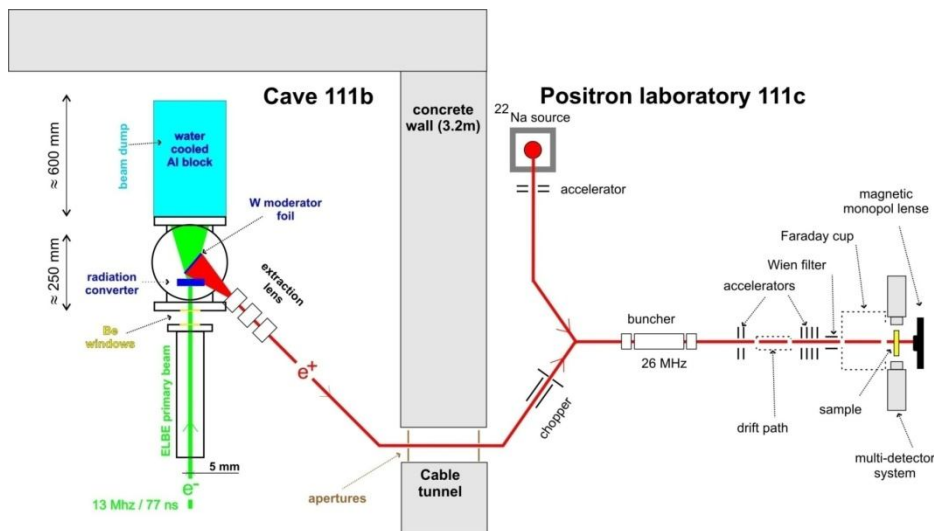

New Concept of EPOS (ELBE Positron Source)

MePS

Monoenergetic Positron Spectroscopy

- Cave 111b / Lab 111d
- monoenergetic (slow) positrons
- pulsed system
- LT, CDBS, AMOC
- Still under construction

Information Depth:
0...5 μm

CoPS

Conventional Positron Spectroscopy

- LT, CDBS, AMOC
- using ^{22}Na foil sources
- He-cryostat
- automated system
- digital detector system

Information Depth:
10...100 μm

GiPS

Gamma-induced Positron Spectroscopy

- Cave 109 (nuclear physics)
- Positron generation by Bremsstrahlung
- Information in complete bulky sample (up to 00 cm^3)
- all relevant positron techniques (LT, CDBS, AMOC)

Information Depth:
0.1 mm ...5 cm

Gamma-induced Positron Spectroscopy



ELSEVIER

Nuclear Instruments and Methods in Physics Research A 495 (2002) 154–160

NUCLEAR
INSTRUMENTS
& METHODS
IN PHYSICS
RESEARCH
Section A

www.elsevier.com/locate/nima

Bremsstrahlung-induced highly penetrating probes for
nondestructive assay and defect analysis


F.A. Selim^{a,*}, D.P. Wells^a, J.F. Harmon^a, J. Kwofie^a, R. Spaulding^a,
G. Erickson^b, T. Roney^c

^aIdaho Accelerator Center, Idaho State University, Campus Box 8263, Pocatello, ID 83209, USA

^bBoise State University, Boise, ID 83725, USA

^cIdaho National Engineering and Environmental Laboratory, Idaho Falls, ID 83415, USA

Received 16 April 2002; received in revised form 13 August 2002; accepted 20 August 2002

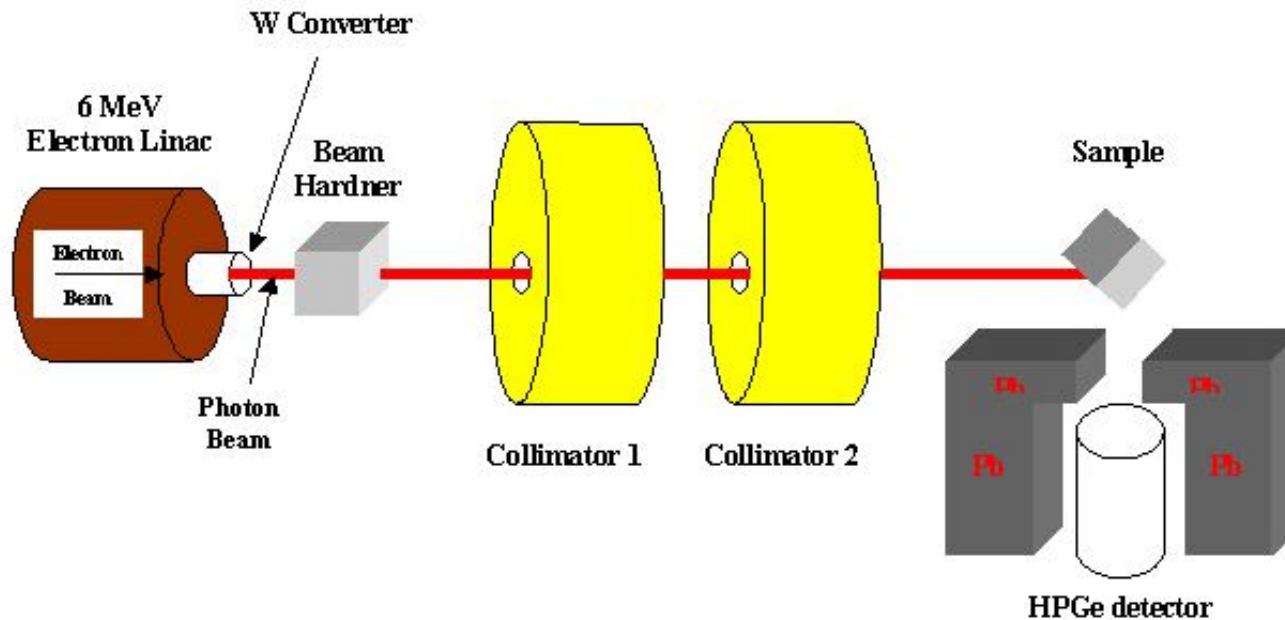

First spectra



Advantages

- information depth 0.1 ... 5 cm; whole sample
- ideal for bulky samples (NDT), liquids, gases, biological objects, coarse powder, dispersions ...

Disadvantage of slow LINACs

- Use of “normal” LINAC with 200 Hz has the problem of high gamma flux in only very few bunches
- Count rate very low, thus no coincidence techniques applicable such as CDBS or AMOC
- Peak / BG ratio bad (1:1)
- no lifetime spectroscopy possible
- Gammas which were scattered into the detector can not be suppressed
- S parameter with huge errors $\pm 4\%$ (more than many effects to be measured)

All this **disadvantages can be overcome** by use of a superconducting LINAC with > 10 MHz

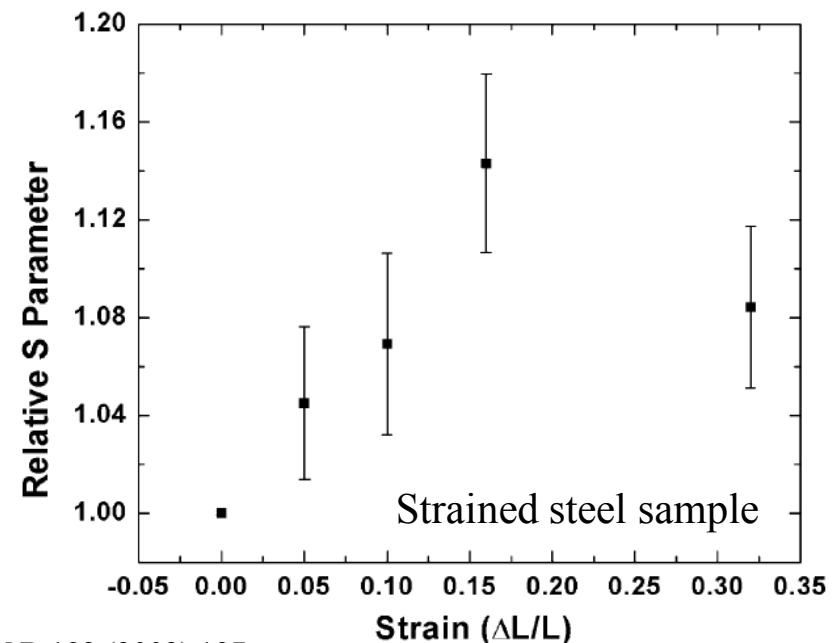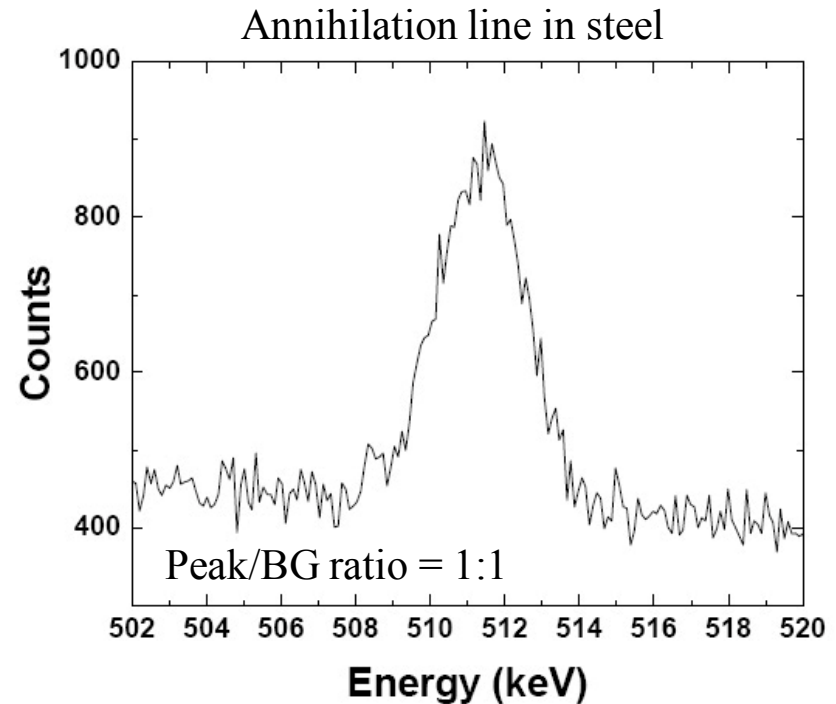

Bremsstrahlung Gamma Source of ELBE (FZ Dresden-Rossendorf)

- Gamma source; pulsed 26 MHz; bunch length < 5 ps; up to 20 MeV; 20 kW beam power
- ideal for GiPS



Available online at www.sciencedirect.com



Nuclear Instruments and Methods in Physics Research A 555 (2005) 211–219

NUCLEAR
INSTRUMENTS
& METHODS
IN PHYSICS
RESEARCH
Section A

www.elsevier.com/locate/nima

The photon-scattering facility at the superconducting electron accelerator ELBE[☆]

R. Schwengner^{a,*}, R. Beyer^a, F. Dönau^a, E. Grosse^{a,b}, A. Hartmann^a, A.R. Junghans^a,
S. Mallion^{a,1}, G. Rusev^a, K.D. Schilling^a, W. Schulze^a, A. Wagner^a

^a*Institut für Kern- und Hadronenphysik, Forschungszentrum Rossendorf, PF 510119, 01314 Dresden, Germany*

^b*Institut für Kern- und Teilchenphysik, Technische Universität Dresden, 01062 Dresden, Germany*

Received 20 July 2005; received in revised form 6 September 2005; accepted 10 September 2005

Available online 10 October 2005

Abstract

A new facility for the production of polarised bremsstrahlung has been built at the superconducting electron accelerator ELBE of the Forschungszentrum Rossendorf. The bremsstrahlung facility and the setup for photon-scattering experiments are designed such that the background radiation caused by the scattering of photons and the production of neutrons is minimised. The sensitive setup in connection with electron energies up to 20 MeV and average currents up to 1 mA delivered by the ELBE accelerator enables novel experiments using photon-induced reactions. First results of photon-scattering experiments are presented.


© 2005 Elsevier B.V. All rights reserved.

PACS: 23.20.-g; 25.20.Dc; 29.30.Kv

Keywords: Superconducting electron accelerator; Polarised bremsstrahlung; Photon scattering


Bremsstrahlung Gamma Source of ELBE (FZ Rossendorf)

- Pulsed gamma source using superconductive Linac ELBE
 - repetition frequency 26 MHz (or smaller by factor 2^n) in CW mode!
 - bunch length < 5 ps
 - up to 20 MeV (we used 16 MeV), no activation of samples by γ -n processes was found
 - average electron current 1 mA = 20 kW beam power; electron beam dump outside lab
 - thus gamma background at target position is very low (Ge detectors with 100% efficiency)
- **Ideal for GiPS ! Is now part of EPOS project – user dedicated positron source.**

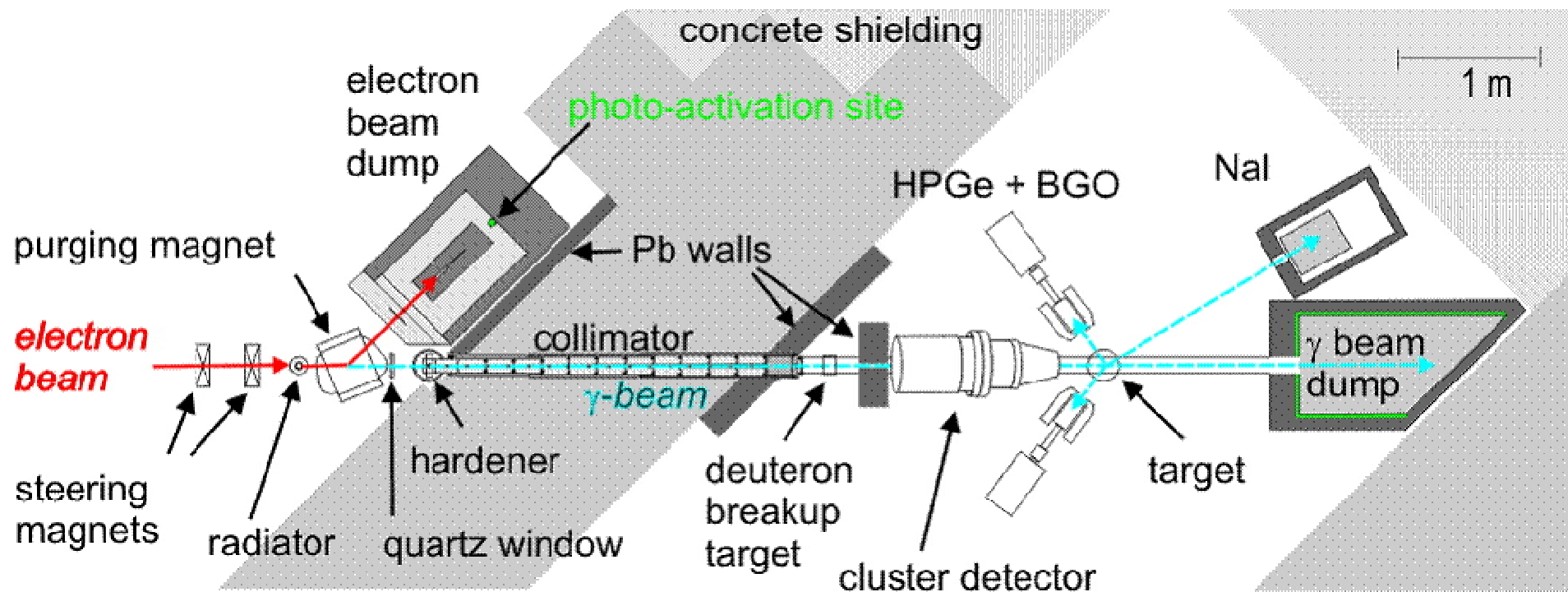

Gamma Spectrum obtained by scattering at ^{208}Pb

- 511 keV annihilation line is by far strongest line although absorber plates were used in front of Ge detectors to reduce the low-energy photons

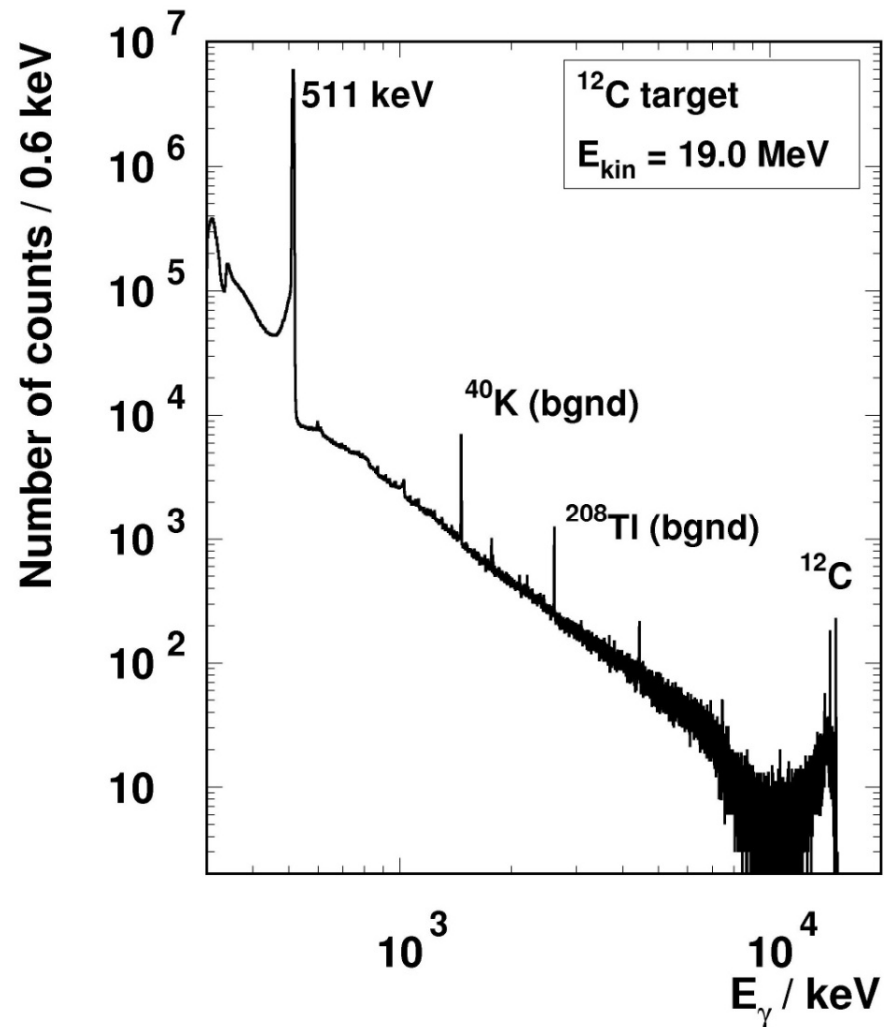


Spectrum of photons scattered from ^{12}C using Bremsstrahlung produced by electrons with an energy of 19 MeV.

Setup extended by BaF₂ detectors for lifetime measurement

- 3 coincident setups were used: 2 AMOC and 1 CDBS spectrometer
- only coincident detection ensures high spectra quality

Problem


- all scattered quanta appear within positron lifetime – time coincidence alone does not reduce background at all
- but distance helps: for 2 x 511 keV quanta in coincidence the distance dependence is proportional to r^{-2}
- for arbitrary scattered gamma it is $\propto r^{-4}$

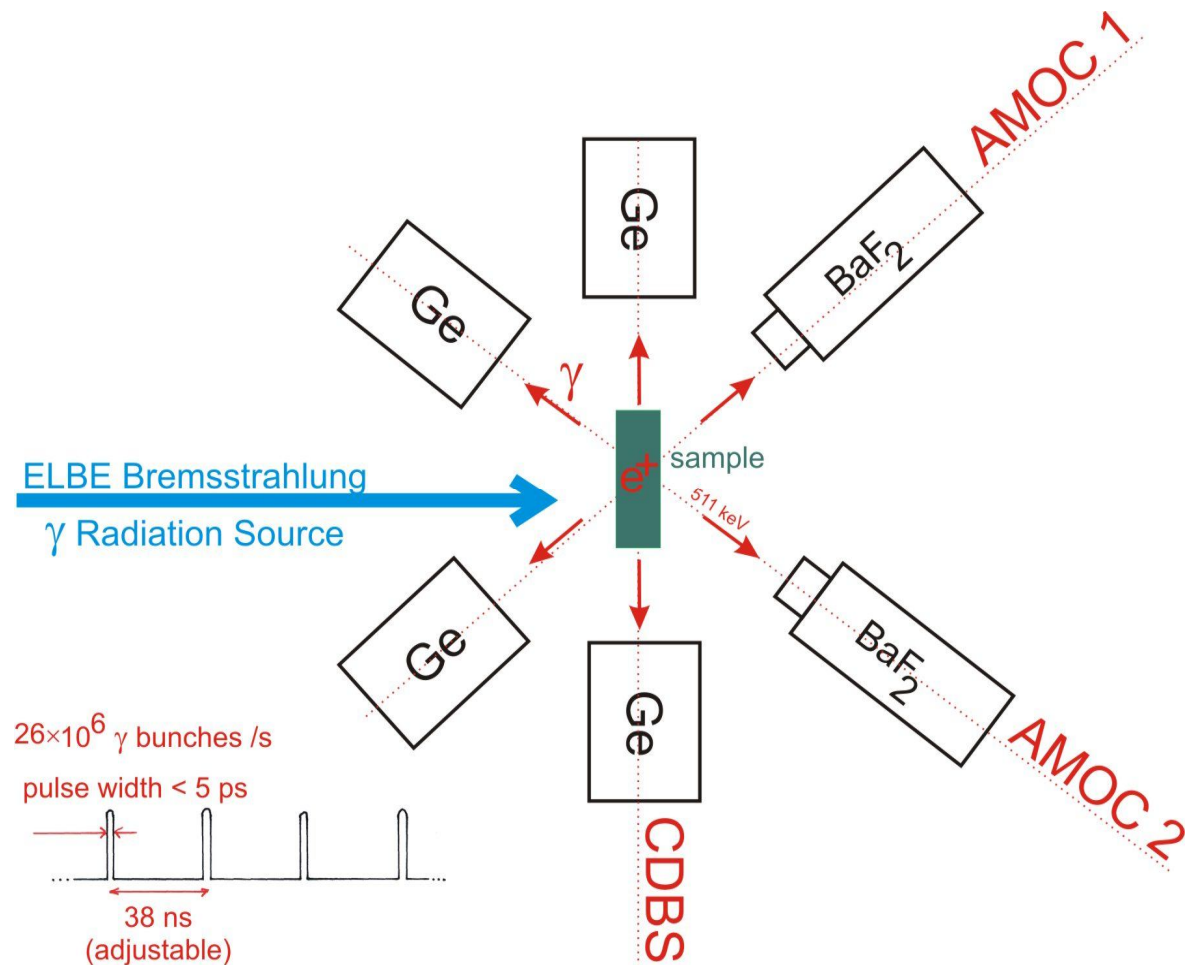


AMOC: Age-Momentum Correlation

CDBS : Concidence Doppler-Broadening Spectroscopy


Setup extended by BaF₂ detectors for lifetime measurement

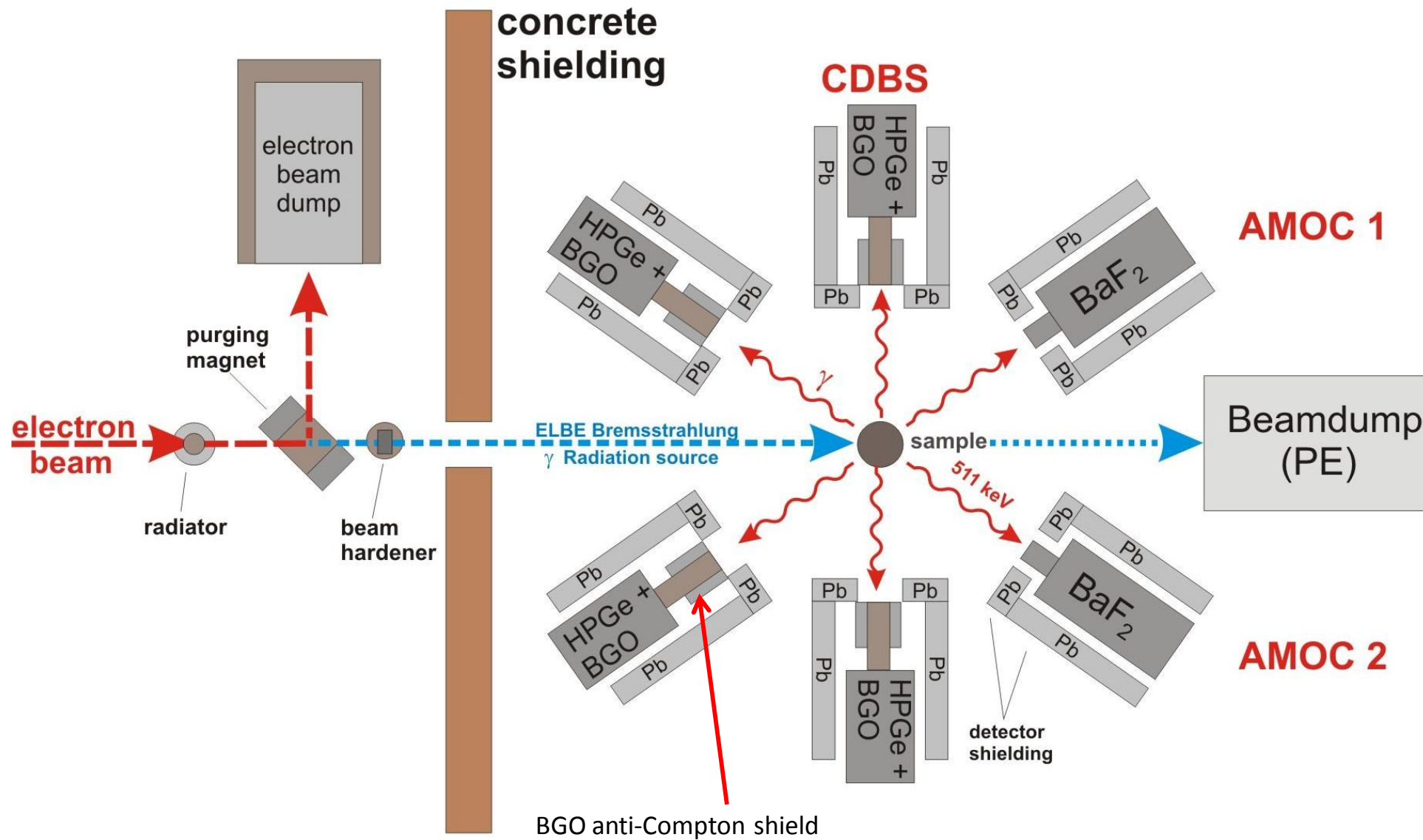

The GiPS setup includes 6 Detectors (4 Ge and 2 BaF₂)

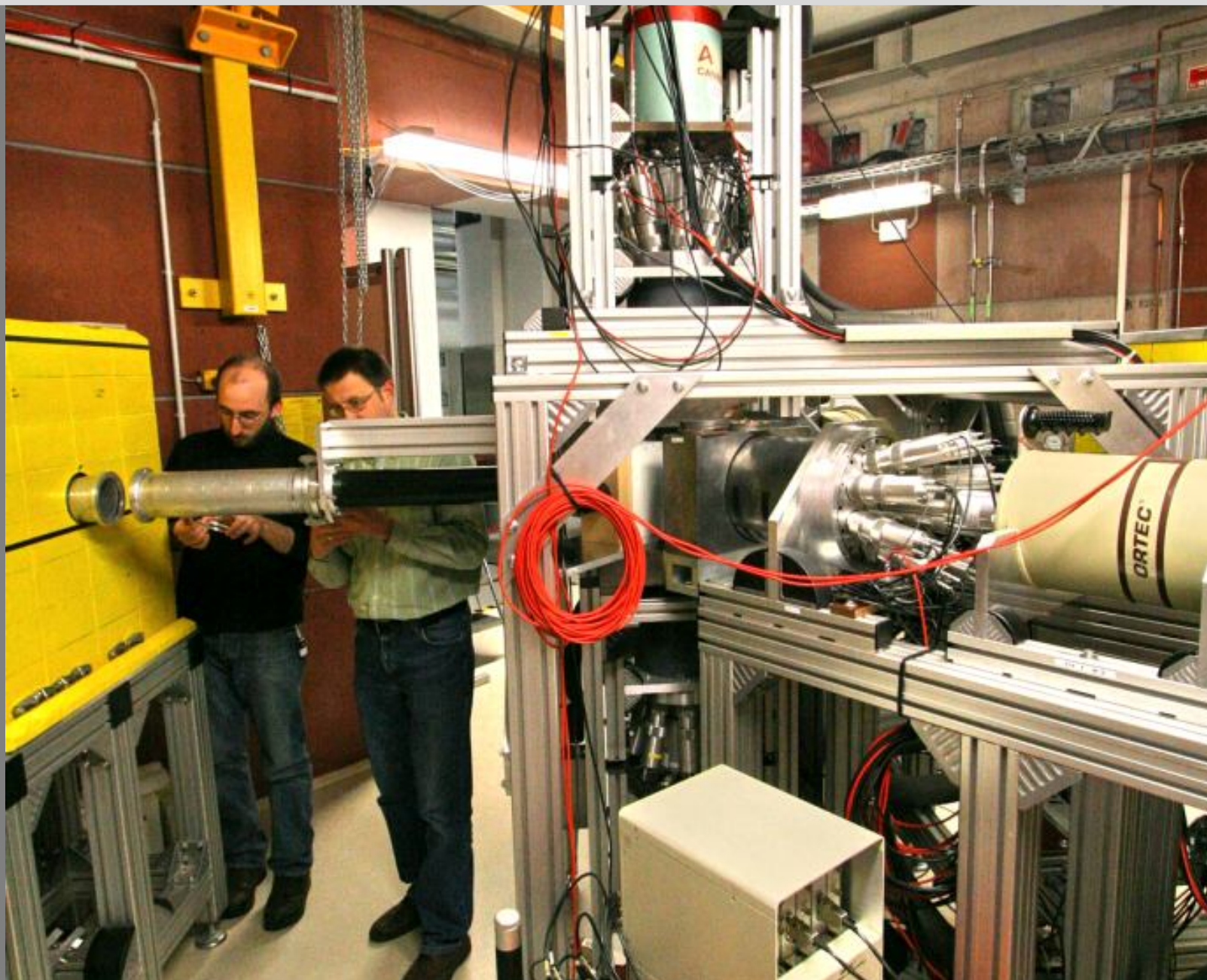


The GiPS setup includes 6 Detectors (4 Ge and 2 BaF₂)

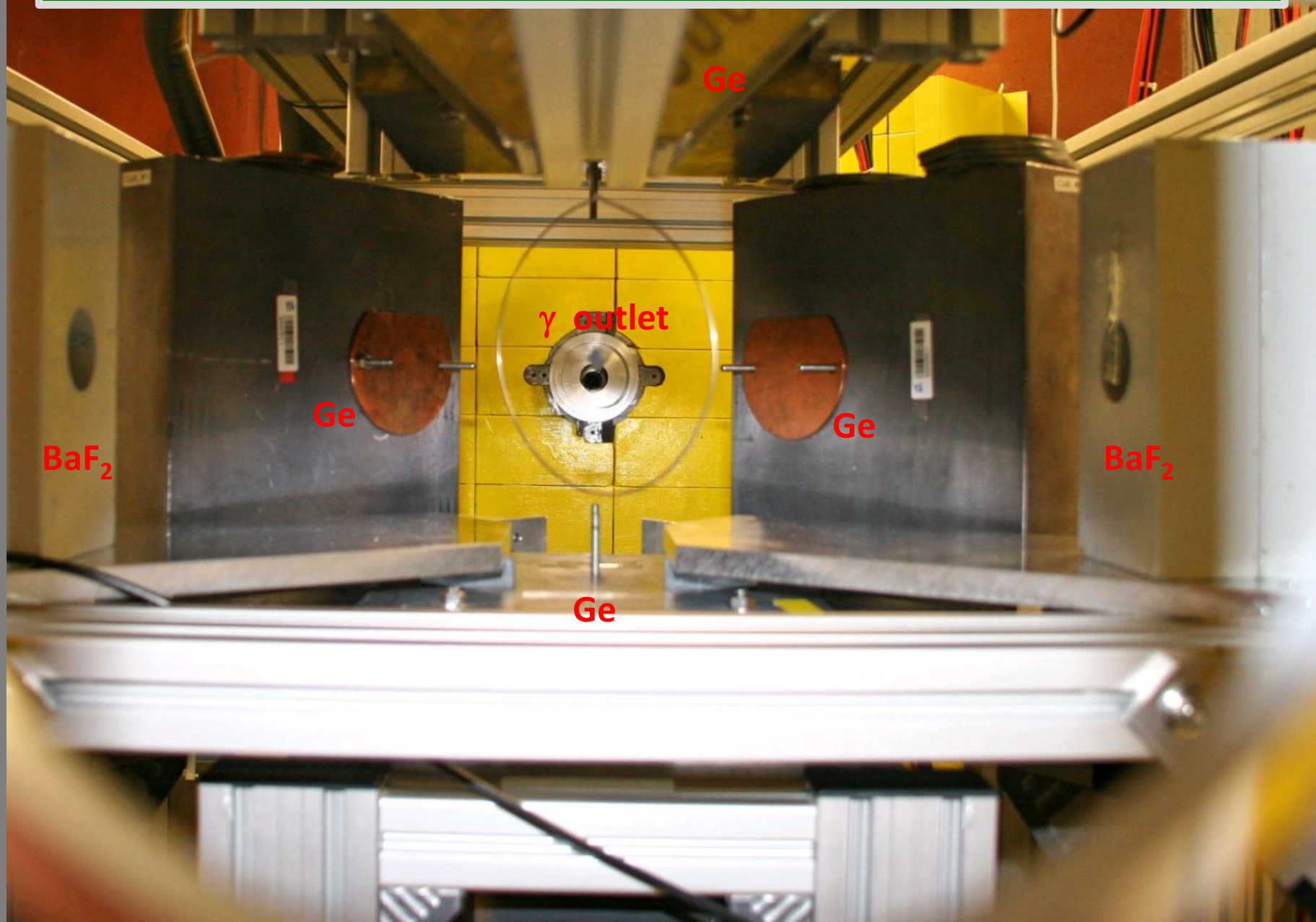


First experiment: annealed and deformed Fe

- Sample 1: well annealed pure sample (< 4N)
- Sample 2: same sample material but 50% thickness reduction by cold rolling

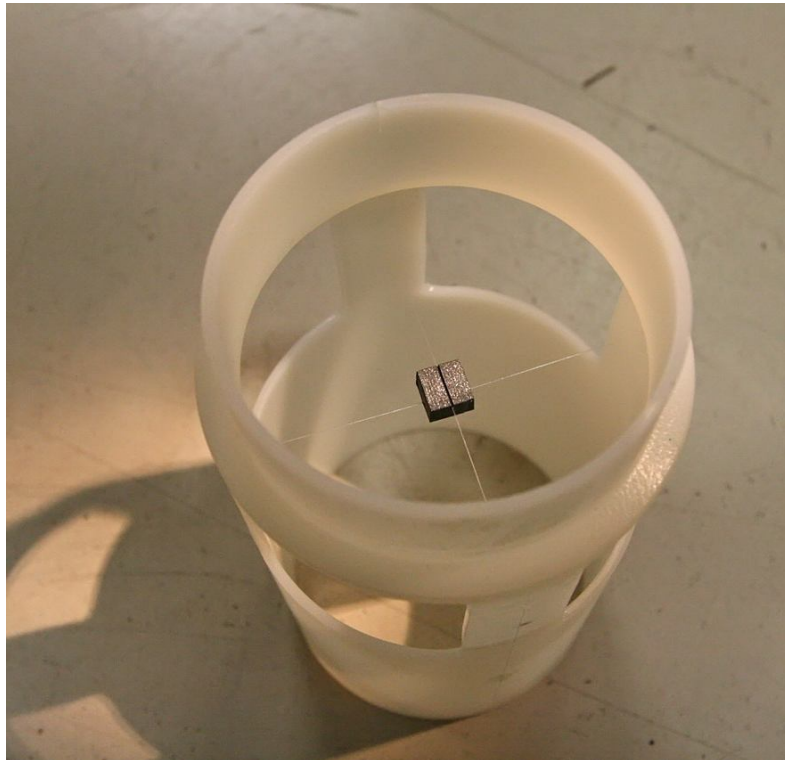

sample holder: Fe sample is hold by thin nylon fibers

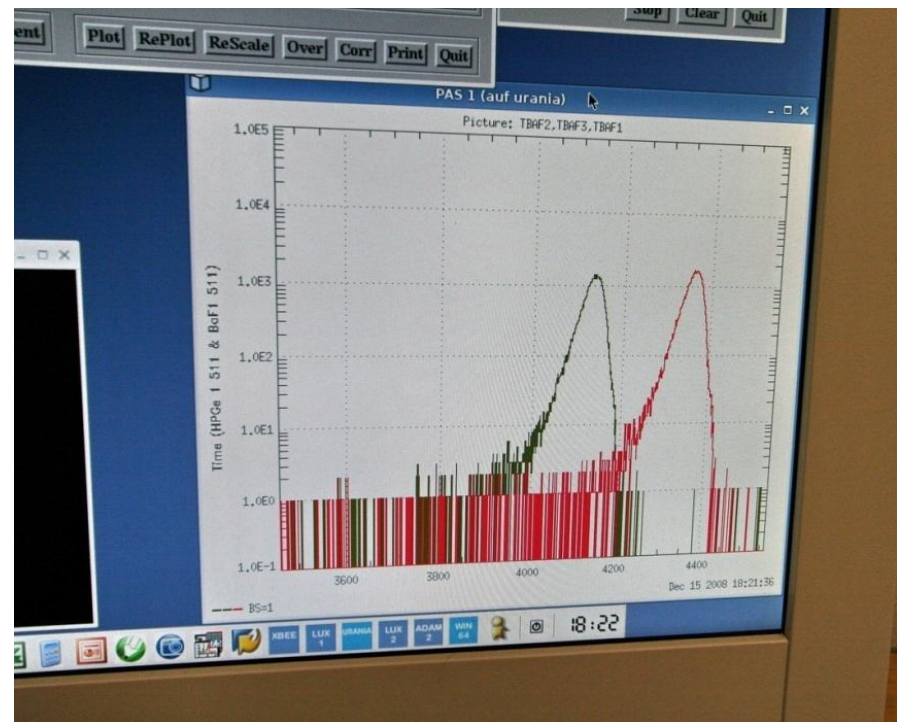

first spectra growing: stop pulse for timing from BaF_2 and start pulse is accelerator bunching pulse of ELBE



Adjustment of single Channel Count Rates

- **all single channel rates together** can be adjusted by changing the electron beam current
- for our Fe samples (0.13 cm^3) we used: $I \approx 100 \mu\text{A}$ – thus even 10 times smaller sample volume would be possible
- here: $8 \times 8 \times 0.5 \text{ mm}^3$ should be enough for having the same rates
- **and the individual rates** can be adjusted:
 - for Ge: by adding absorber materials (here 5 mm Cu plate)
 - for BaF₂: increase of distance
- Read the **counter display** as:
 - first 4 rates from top: 4 Ge detectors (1 and 3 for the AMOC system, 2 and 4 for CDBS)
 - rate 5 and 6: rates of BaF₂ detectors
 - rate 7 and 8: LT coincidence rates of both AMOC systems
- Thus a single channel ...
 - ...lifetime spectrum (4×10^6 counts) needs about 7 min, there were 2 detectors, thus 3.5 min
 - ...Doppler spectrum (1×10^6 counts) requires 90 s (50% of whole intensity is within 511 peak) – but there are 4 Ge detectors, thus 20 s should be enough




Single-channel Ge Spectrum of annealed Fe

- count rate about 20 kHz (200 kHz would be theoretically possible); total counts in example: 8×10^6
- about 50% of intensity in 511 peak
- decrease below 350 keV due to 5 mm Cu absorber plates in front of Ge detectors
- detection with analog electronics

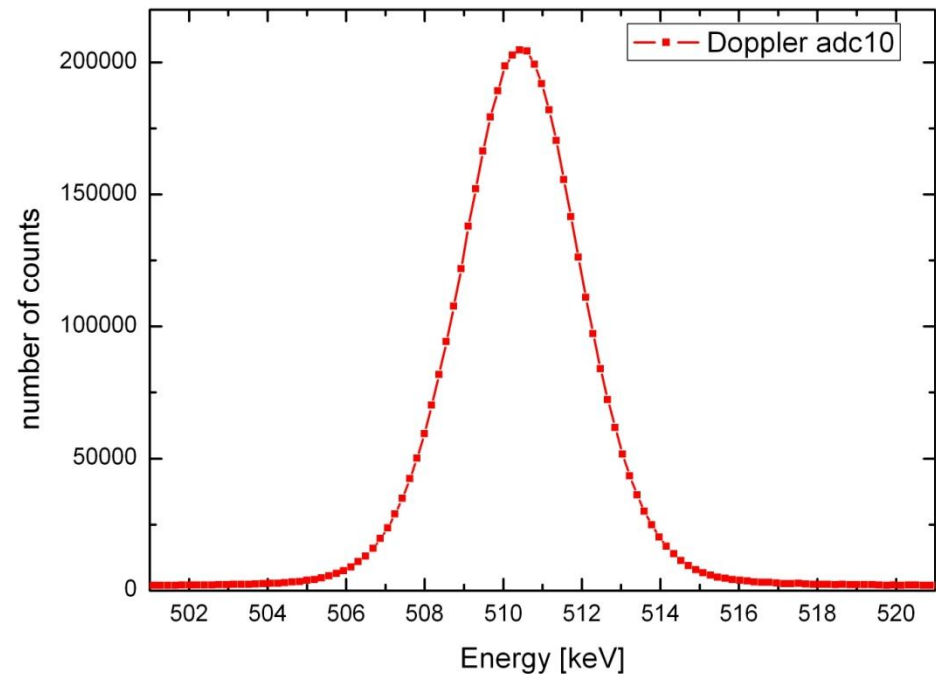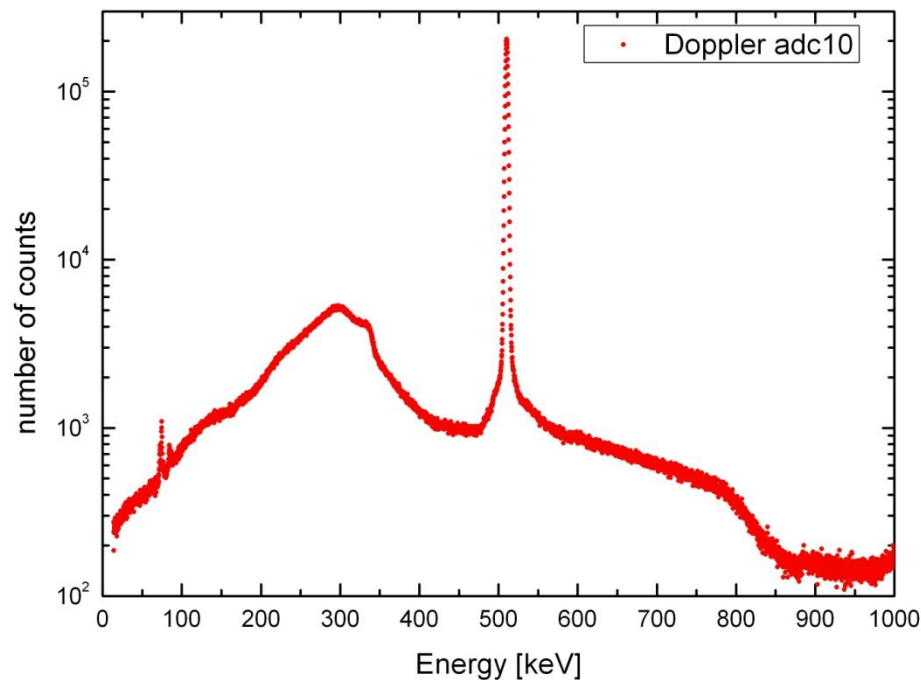

Comparison annealed and deformed Fe

- expected behavior
- curve of deformed Fe is distinctly taller due to open-volume defects and thus increased fraction of annihilation with valence electrons (small energies – small Doppler shift)

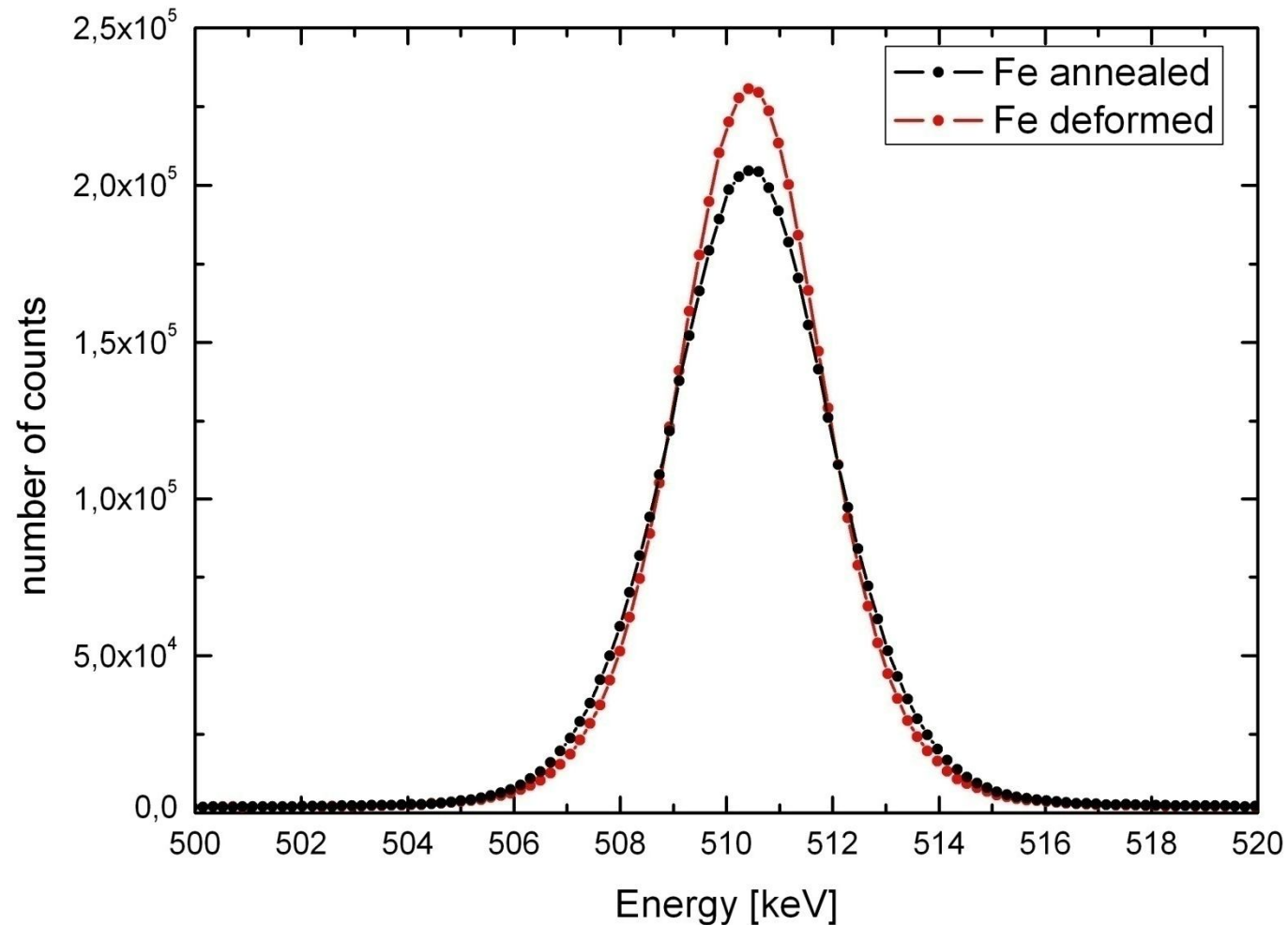

Single-channel Lifetime Spectrum of Fe

- single channel lifetime spectrum (Start: BaF₂ - Stop: accelerator machine pulse)
- obvious: without coincidence – many quanta different from annihilation quanta are detected
- scattered photons
- they are not visible in Doppler spectrum - appear in normal 511 keV peak
- only coincidence measurement ensures suppression of scattered gamma (next slide)

delay of side peaks



1: -0.84 ns = 25.2 cm

2: -0.42 ns = 12.6 cm

3: 1.8 ns = 54 cm


4: 7.0 ns = 2.1 m

scattering occurs from all over the detectors, even from the γ beam dump

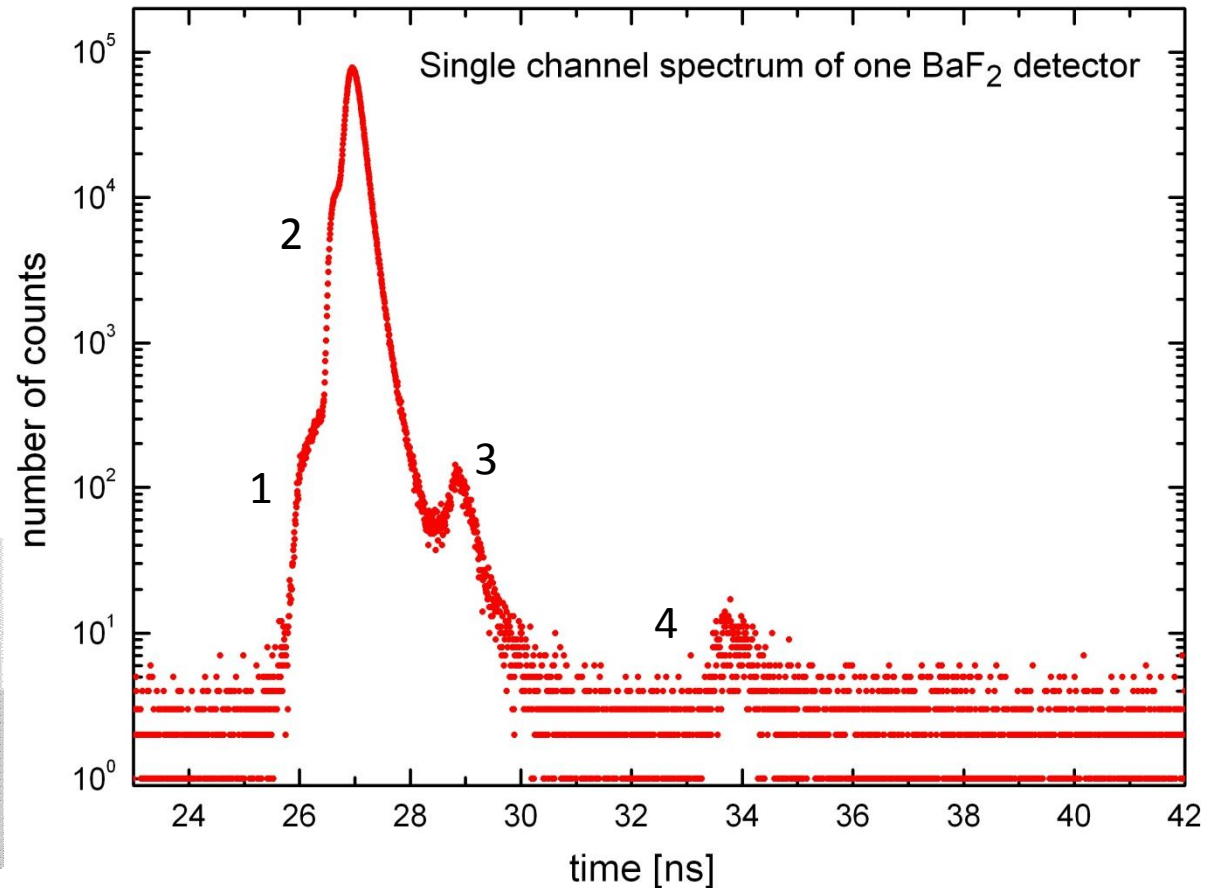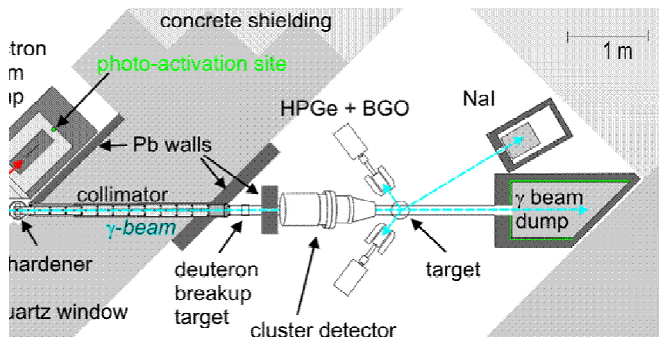


Coincident lifetime spectrum: annealed Fe

- here coincidence with Ge detector
- spectrum is projection to the time scale of AMOC spectrum
- Count rate for single AMOC spectrum = 160 /s
- both spectra: 320 /s – 1 million counts in 1 hour
- Time resolution = 210 ps
- BG/Peak = 1.7×10^{-5}
- 350 ps & 1.5 ns: annihilation at vacuum tube (polyethylene)
 - two areas below t_0 exhibit sometimes disturbances
 - probably due to detection of scattered gammas (at sample and vacuum tube)
 - but nothing after t_0

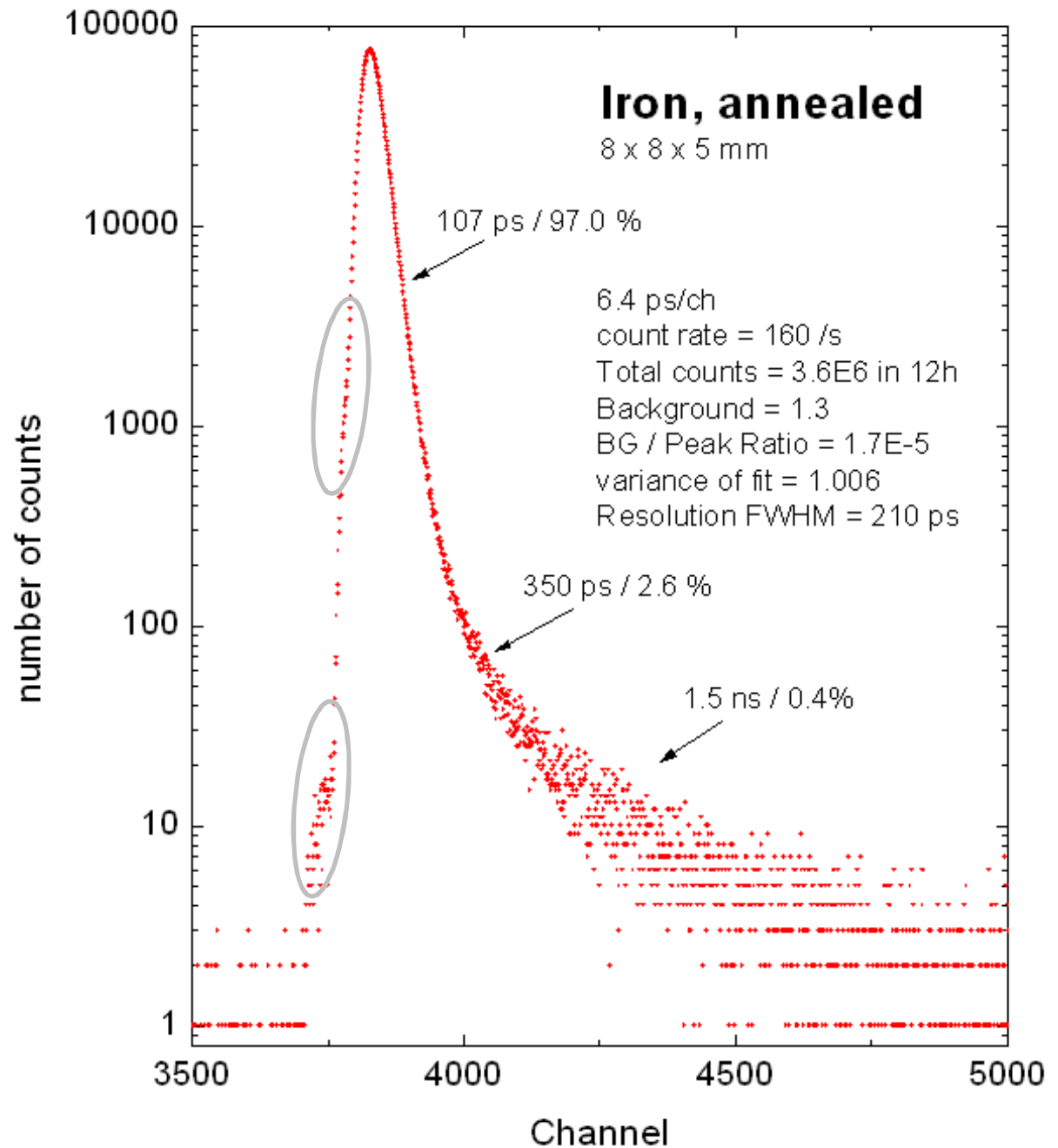


Residuals of fit show perfect fit

- analysis by LT 9.0 (J. Kansy)

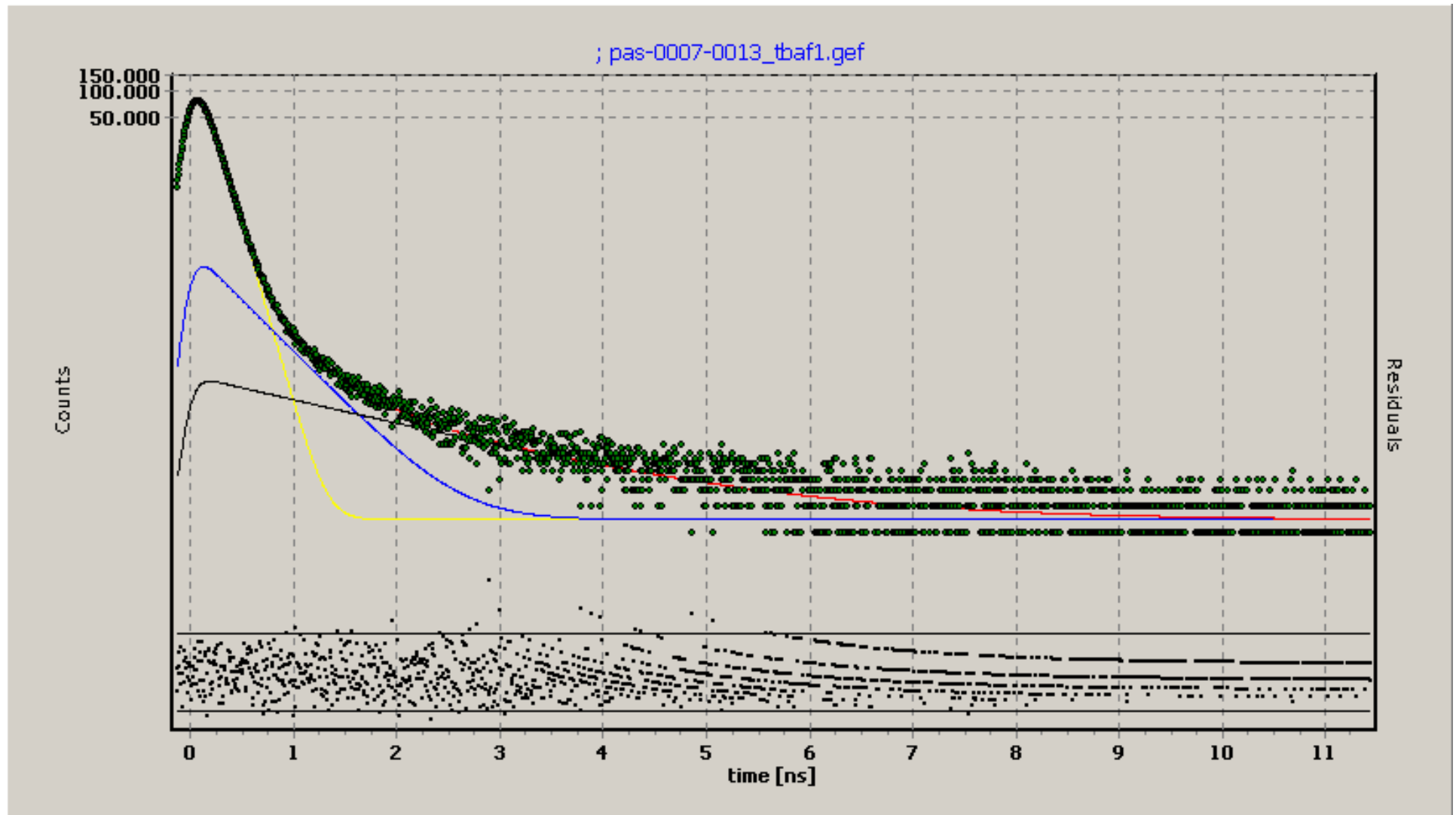


Comparison: GiPS spectrum with conventional measurement

- same sample material – almost same statistics, similar time resolution
- conventional measurement with ^{22}Na source 20 μCi (0.7 MBq) in sandwich geometry
- advantage of periodic positron source is obvious: background distinctly reduced
- result of spectra analysis is the same: 107 ps (bulk value for Fe; corresponds to literature)

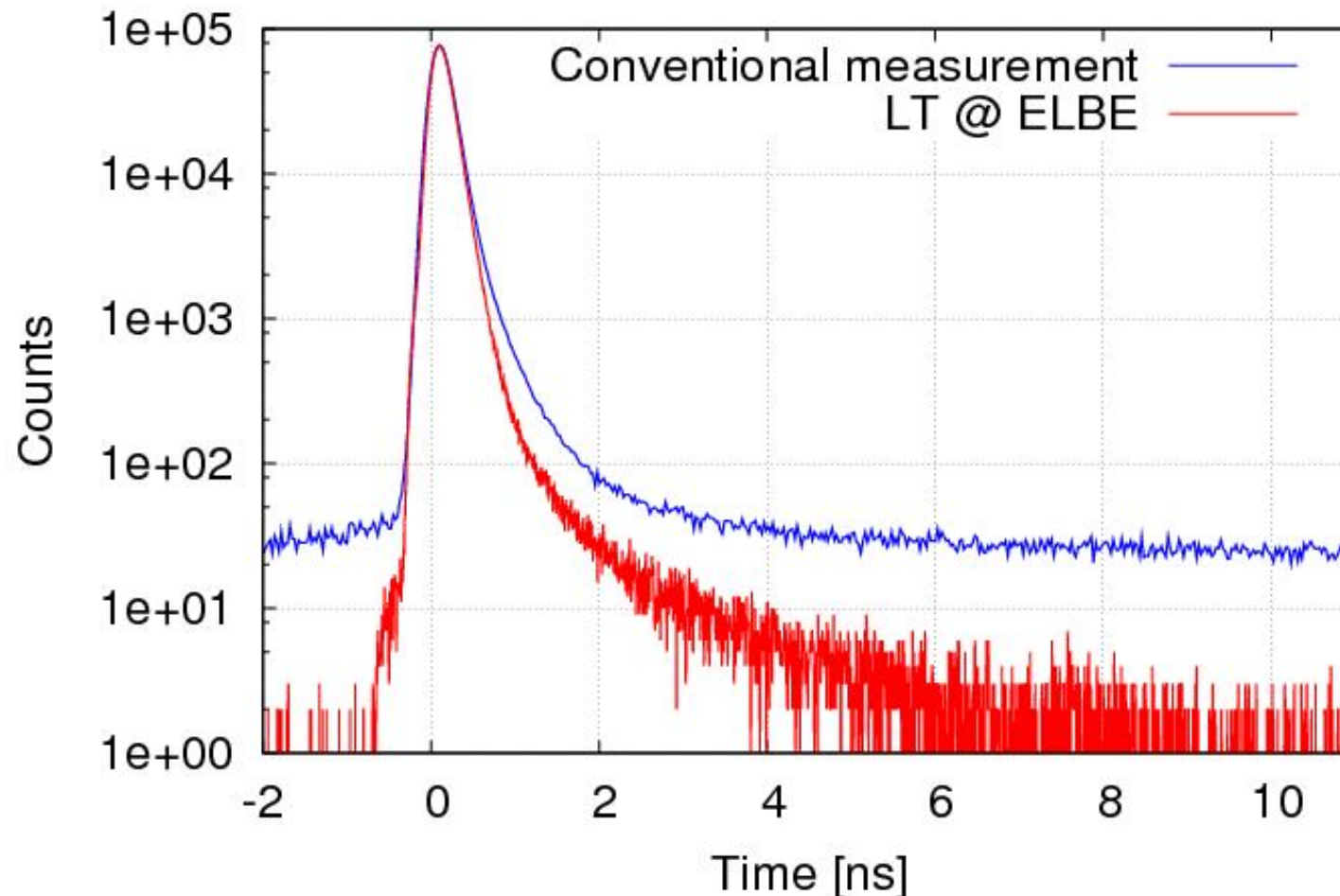

Comparison annealed and deformed Fe


- two mechanically identical samples were prepared
- Fe annealed (1100°C; 2h in vacuum) and Fe (50% thickness reduction by cold rolling)
- spectra were easily decomposed
- expected results: annealed sample – one component 107 ps; deformed sample has 158 and 401 ps (dislocations and small vacancy clusters)

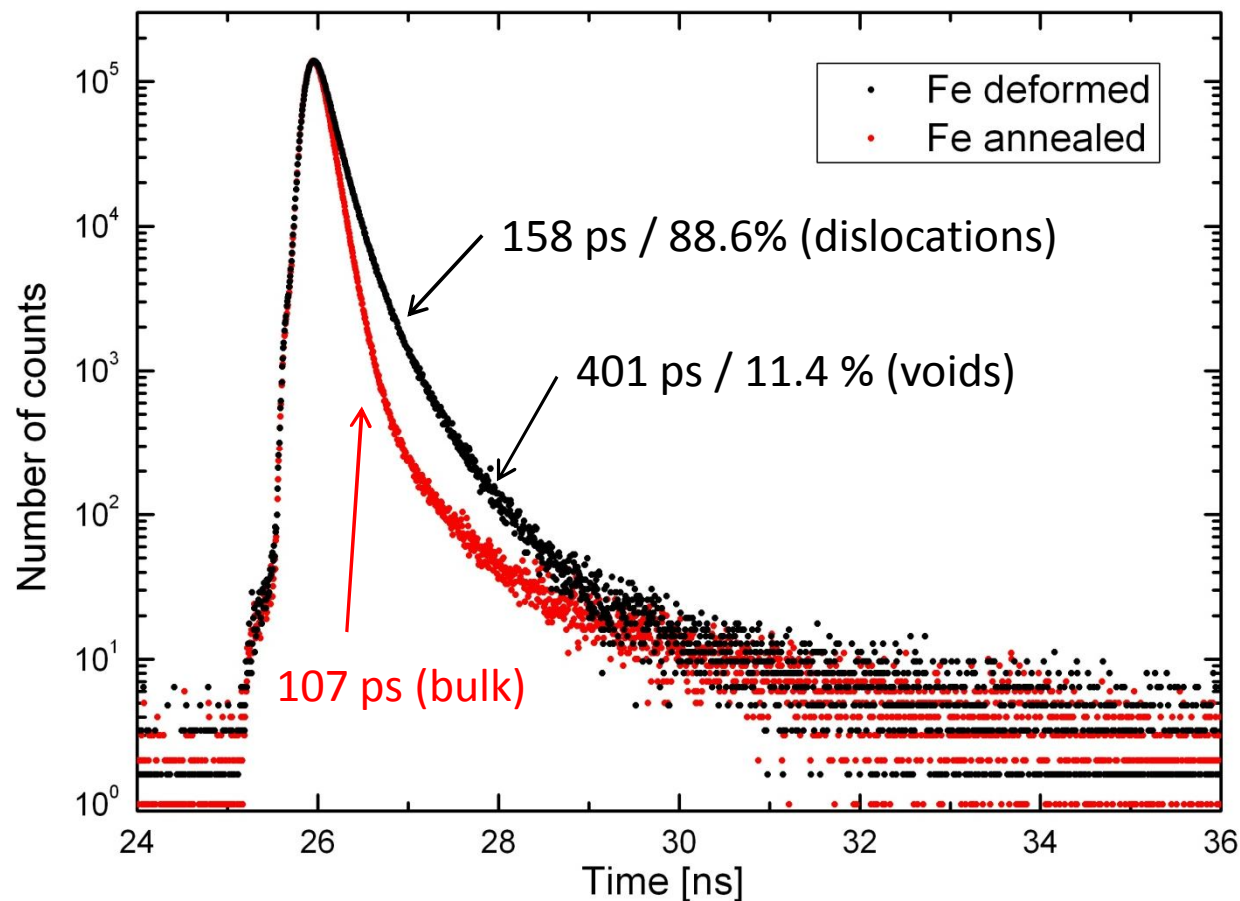


AMOC spectrum of annealed Fe

- AMOC: measurement of momentum of annihilating electron as function of positron age
- AMOC detection is not an extra gimmick, but is required to maintain quality of spectra
- only by coincident measurement of 511 keV annihilation line: suppression of scattered gamma (can be concluded from lifetime spectra)

2D-Plot

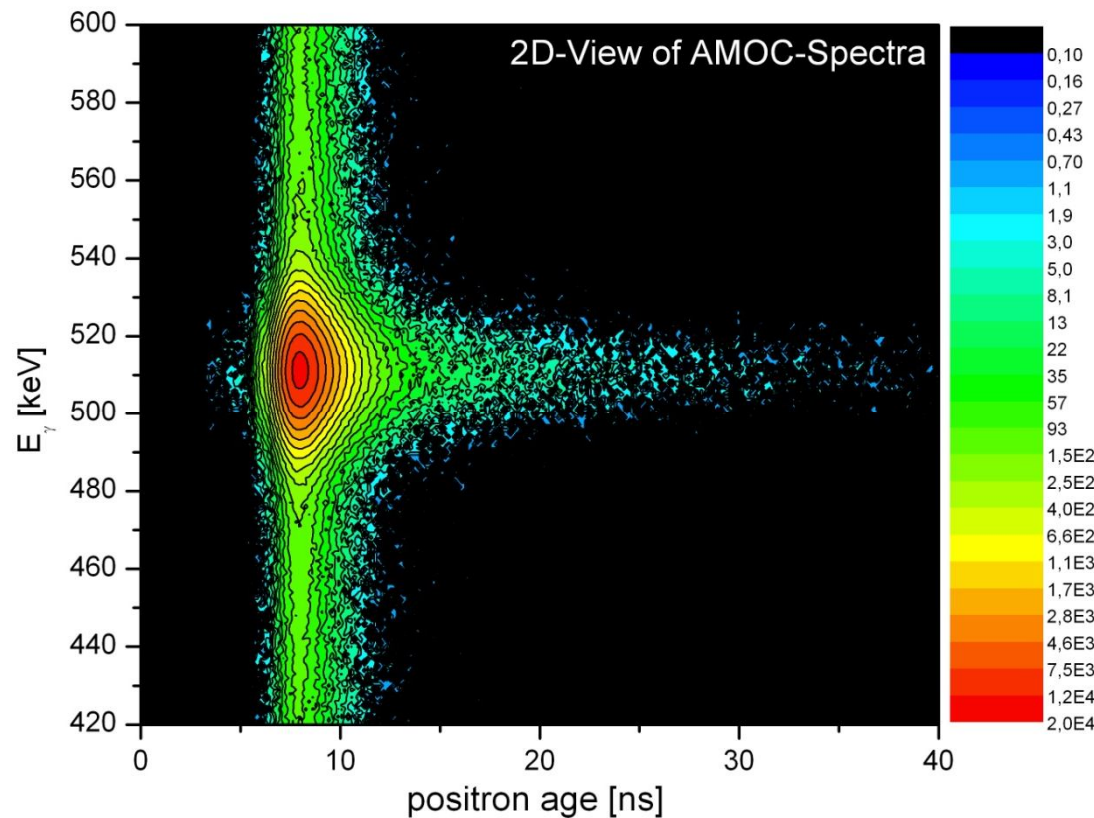


3D-Plot

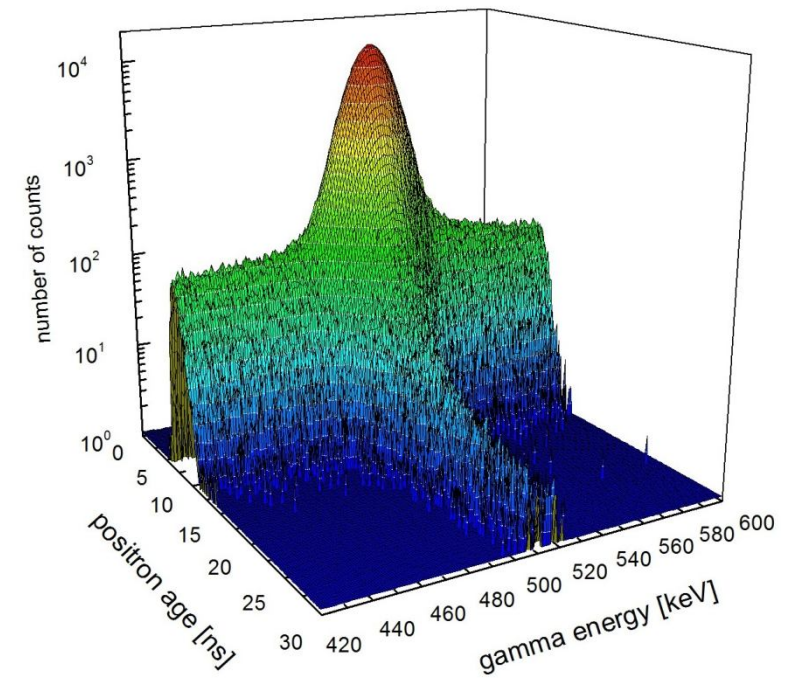



Coincidence Doppler-Broadening Spectroscopy of Fe sample

2D-View of CDB-Spectra

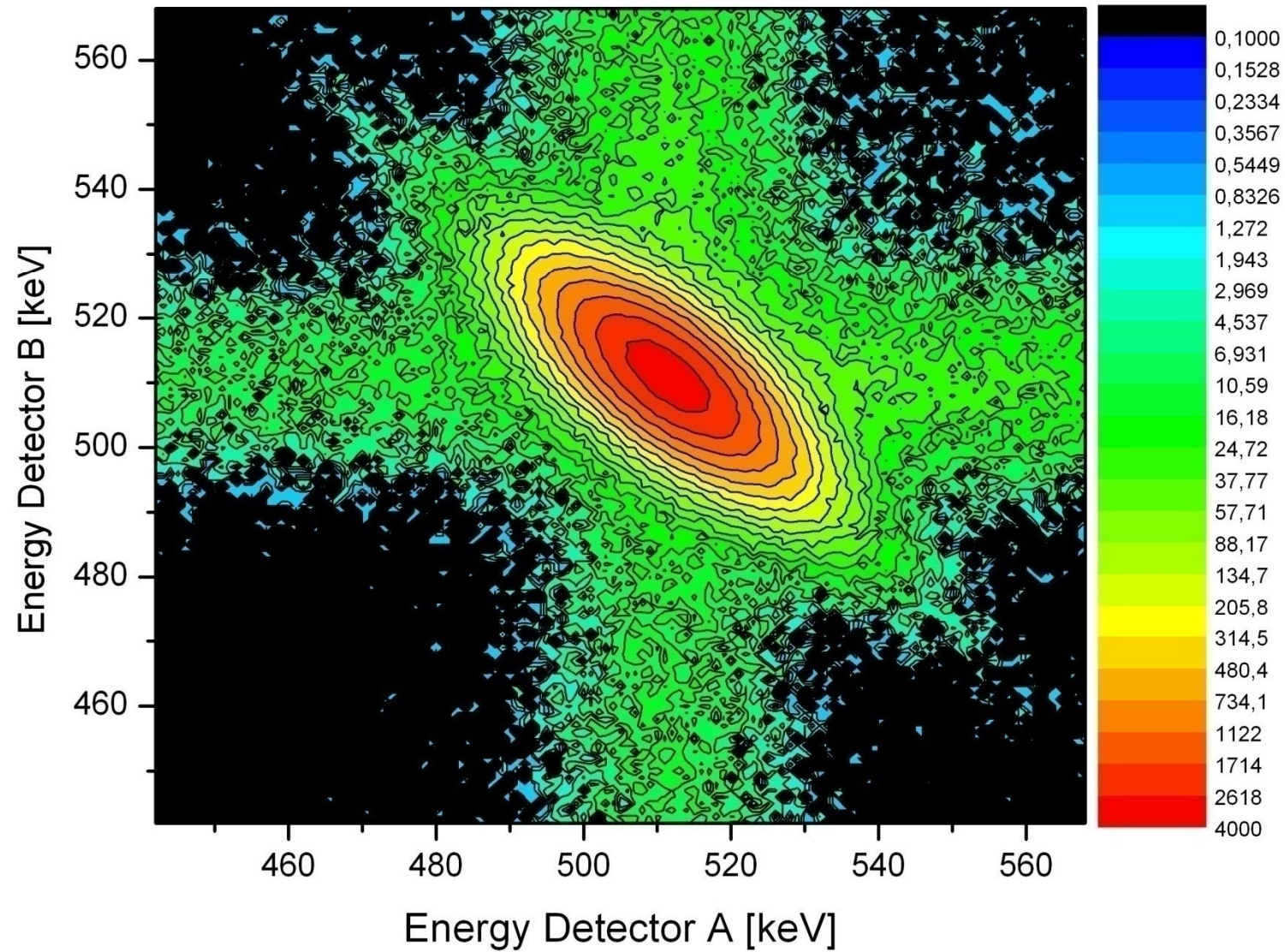


Sample holder for liquids

- Water is hold between two thin Kapton foils (few μm)
- inner diameter of Al rings is 50 mm
- Kapton was used because of simple lifetime spectrum (single component of 382 ps)
- but volume fraction of Kapton is very small (0.4%): almost no contribution to spectrum expected
- gamma beam has 20 mm diameter

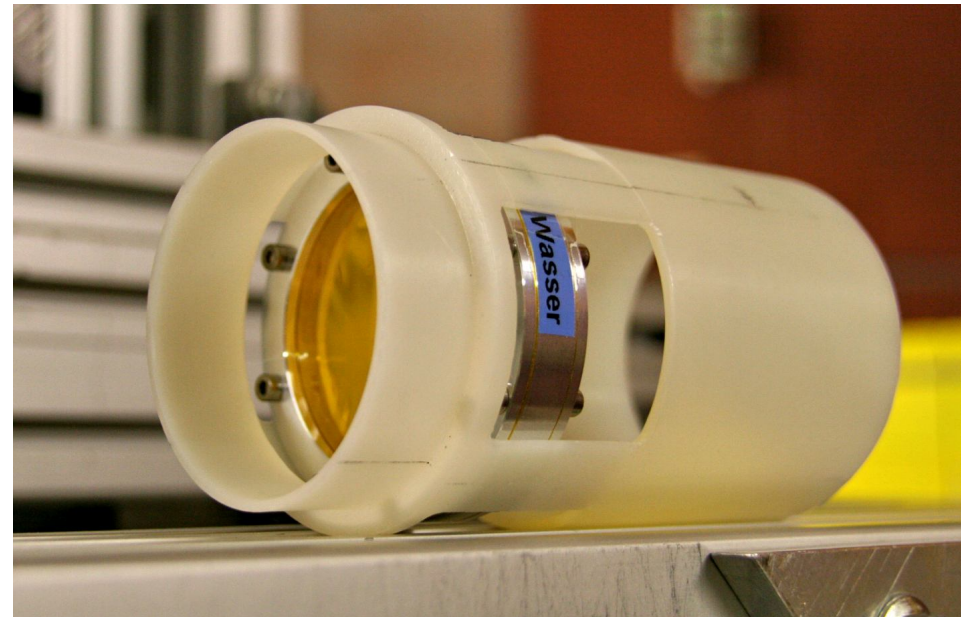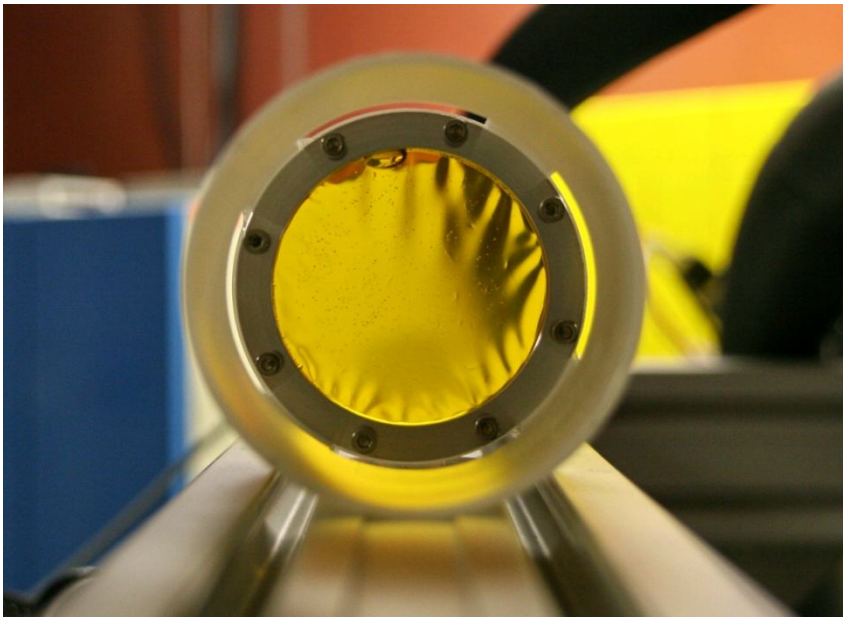


Water at RT

- total count rate: 5×10^5
- no such visible deviations on $t < t_0$ like for Fe (due to much smaller gamma scattering compared to Fe)

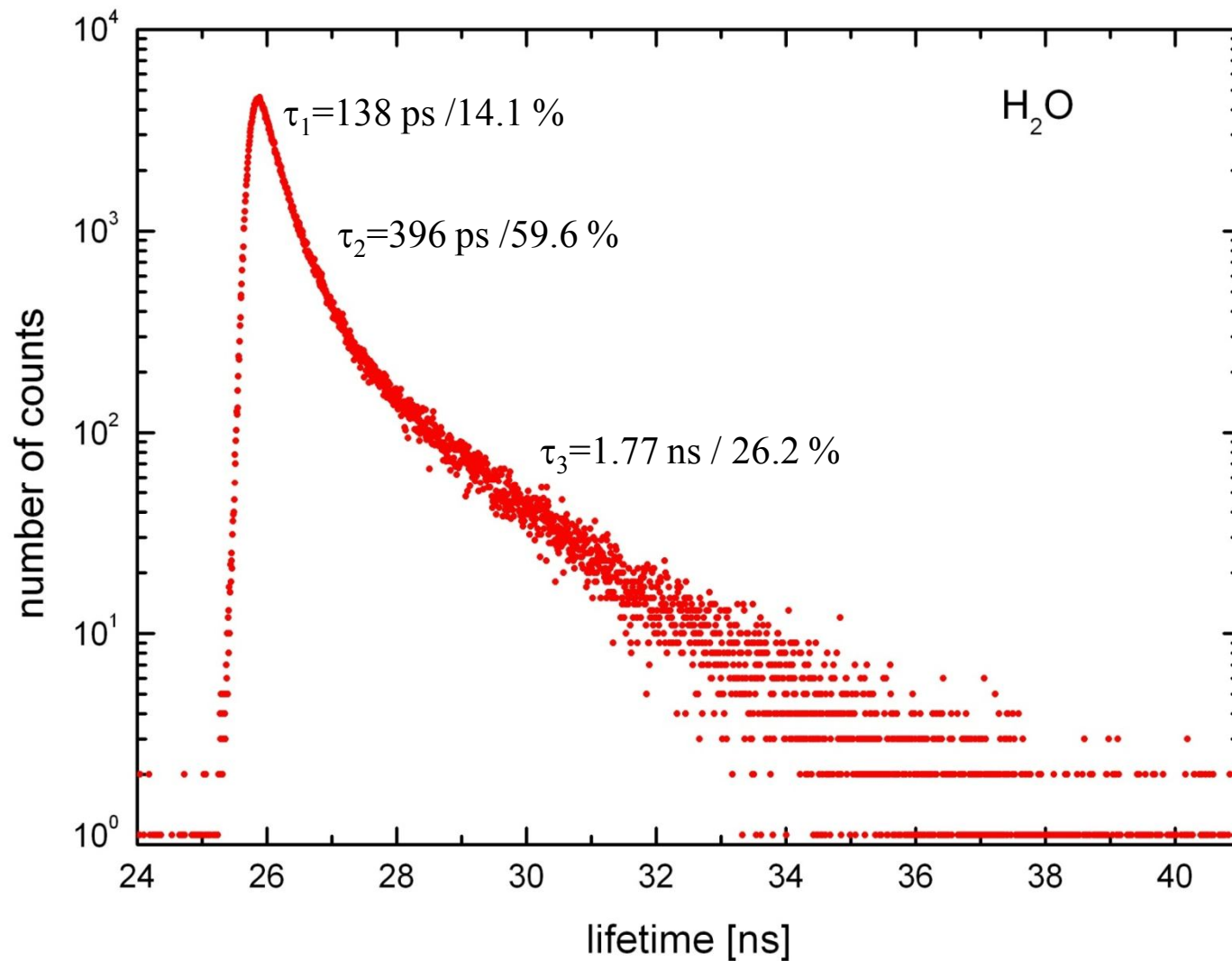

Comparison of Water Spectrum with conventional setup

- total count rate: $12 \times 10^6 \text{ s}^{-1}$

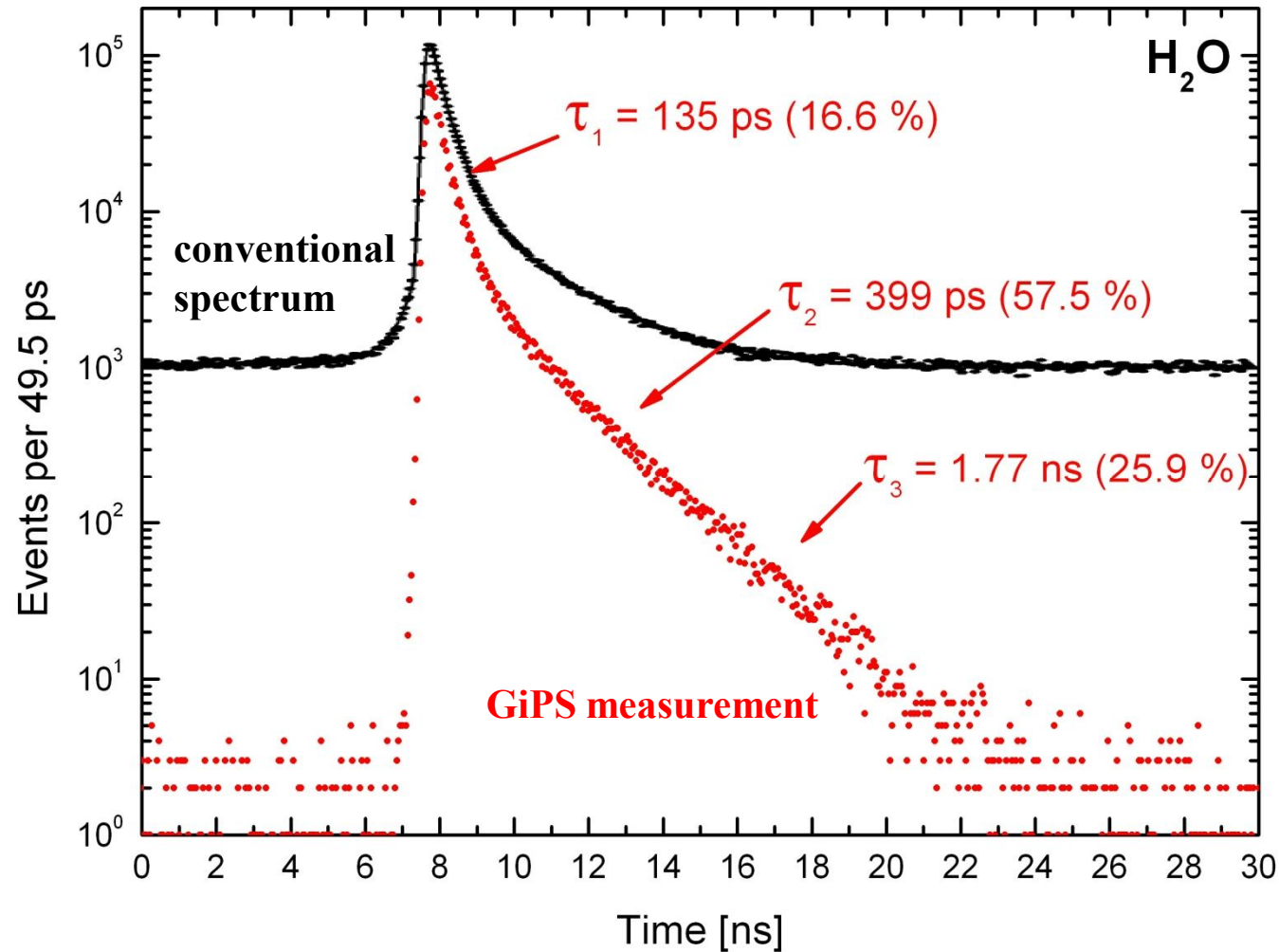





- Black spectrum: conventional measurement by K. Kotera, T. Saito, T. Yamanaka, Phys. Lett. A 345, (2005) 184

Amorphous Silica Glass

- round piece 1.5 cm thick, about 5 cm³
- lifetime spectrum: total count rate: 2x10⁶
- same sample was measured conventionally in 1978 also in the same institute (former ZfK Rossendorf):

151 ps - 523 ps - 1.57 ns (FWHM \approx 350 ps)

G. Brauer et al., Appl. Phys. 16 (1978) 231

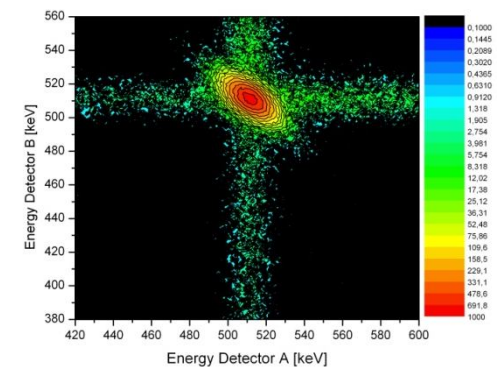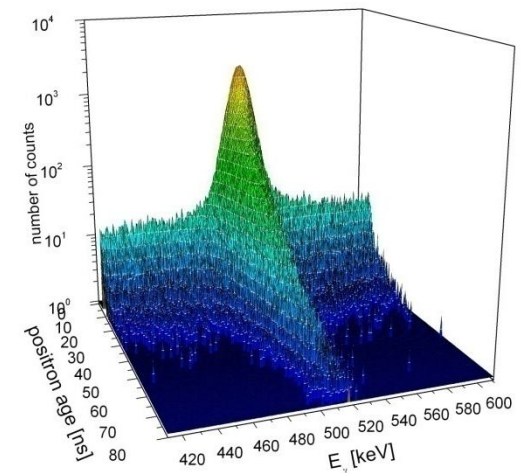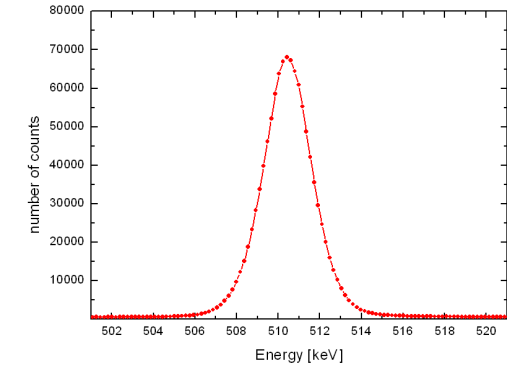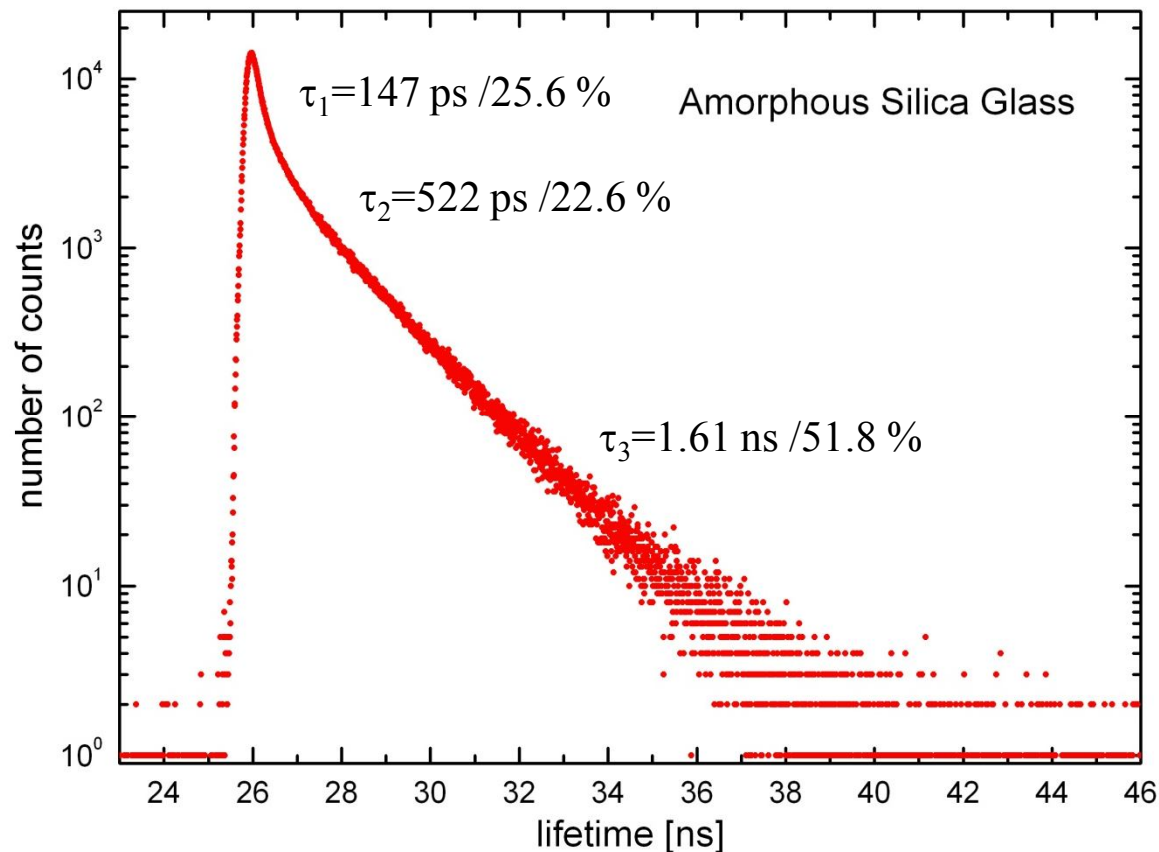

Conclusions

- new concept of EPOS project is now extended to use mono-energetic Positrons (MePS), Gamma-induced (GiPS) and conventional spectroscopy (CoPS)
- all spectrometers are equipped with LT, CDB, AMOC
- fully digital system (in the future)
- EPOS can cover sample thickness range from 10 nm to 10 cm (7 orders of magnitude)
- MePS still under construction
- GiPS has been tested successfully
 - GiPS only possible because of the unique properties of the ELBE Linac (cw mode of 26 MHz intense and extremely short electron bunches, < 5ps bunch length)
 - background suppression by coincident measurement of Lifetime and Doppler (AMOC)
 - surprisingly good spectra quality
 - coincidence between 2 BaF₂: resolution improves by 24% (FWHM = 160 ps)
 - problem: heating / cooling of sample because in holder positrons are also generated

Talk available at <http://positron.physik.uni-halle.de>